

Operating instructions

Refrigerator/freezer

LABO series



 **Kirsch**
Made in Germany. Since 1865.

Read instruction manual before starting work!

Philipp KIRSCH GmbH
Im Lossenfeld 14
77731 Willstätt-Sand
GERMANY
Telephone: +49 781 9227-0
Fax: +49 781 9227-200
Email: info@KIRSCH-medical.de
Internet: www.KIRSCH-medical.de

D1006_05.21_Operating instructions LABO version 2, en_GB

About this instruction manual

This instruction manual was created for the product “Refrigerator/freezer” (hereafter referred to as “unit”).

Persons who work with the unit must have carefully read and understood this instruction manual before any work begins. To ensure safe working conditions, all specified safety warnings and instructions in this instruction manual must be observed.

In addition, special owner obligations may apply.

Keeping the manual

This instruction manual makes it possible to handle the unit safely and efficiently. This instruction manual is part of the unit; it must be kept in the immediate vicinity of the unit and be accessible to staff at all times.

Target audiences

This instruction manual is designed to provide information to the following target audiences:

- Owner of the unit
- Users of the unit

A separate service manual for this unit is available for technical customer service (hereafter referred to as “service department”).

Models

This manual is valid for the following units:

Model	Factory number from
LABO 100 PRO-ACTIVE	100 05 25000 / 100 33 25000
LABO 126 PRO-ACTIVE	125 05 25000 / 125 34 25000
LABO 288 PRO-ACTIVE	280 15 25000 / 280 34 25000
LABO 340 PRO-ACTIVE	330 12 25000 / 330 31 25000
LABO 468 PRO-ACTIVE	460 05 25000 / 460 31 25000
LABO 520 PRO-ACTIVE	500 01 25000 / 500 30 25000
LABO 720 PRO-ACTIVE	700 01 25000 / 700 30 25000
LABO 520 CR PRO-ACTIVE	500 02 25000 / 500 42 25000
LABO 720 CR PRO-ACTIVE	700 02 25000 / 700 36 25000
LABO 720 CHROMAT PRO-ACTIVE	700 20 25000 / 700 32 25000
LABO 288 ULTIMATE	280 70 25000 / 280 80 25000
LABO 340 ULTIMATE	340 70 25000 / 340 80 25000
LABO 468 ULTIMATE	460 70 25000 / 460 80 25000
LABO 520 ULTIMATE	500 70 25000 / 500 80 25000
LABO 720 ULTIMATE	700 70 25000 / 700 80 25000

Model	Factory number from
LABO 720 CHROMAT ULTIMATE	700 74 25000 / 700 87 25000
FROSTER LABO 330 PRO-ACTIVE	320 01 25000
FROSTER LABO 530 PRO-ACTIVE	500 08 25000
FROSTER LABO 730 PRO-ACTIVE	700 04 25000
FROSTER LABO 530 CR PRO-ACTIVE	500 09 25000
FROSTER LABO 730 CR PRO-ACTIVE	700 05 25000
FROSTER LABO 330 ULTIMATE	320 80 25000
FROSTER LABO 530 ULTIMATE	500 84 25000
FROSTER LABO 730 ULTIMATE	700 84 25000

Illustrations

Illustrations in this manual are designed as an aid to basic comprehension and may deviate from the version at hand.

Manufacturer's address

Manufacturer	Philipp Kirsch GmbH
Address	Im Lossenfeld 14
	77731 Willstätt-Sand
	GERMANY
Telephone	+49 781 9227-0
Fax	+49 781 9227-200
Email	info@kirsch-medical.de
Internet	www.kirsch-medical.de

Business hours:

- Monday to Thursday: 8:00 am to 12:15 pm, 1:15 pm to 4:30 pm
- Friday: 8:00 am to 12:00 pm, 1:00 pm to 4:00 pm

Service contact

Service department	Philipp Kirsch GmbH
Address	Im Lossenfeld 14
	77731 Willstätt-Sand
	GERMANY
Telephone	+49 781 9227-777
Fax	+49 781 9227-200
Email	info@kirsch-medical.de
Internet	www.kirsch-medical.de

Orders are accepted during business hours.

More information

If you have questions or comments regarding this instruction manual or the unit, please contact your authorised regional specialist dealer or contact KIRSCH directly.

Table of contents

1	Product description	8
	1.1 Unit overview.....	8
	1.2 Display and control elements.....	12
	1.2.1 Design of the display and control unit.....	12
	1.2.2 Function of the key switch	12
	1.2.3 Function of buttons and displays.....	12
	1.3 Scope of delivery.....	14
	1.4 Interfaces.....	14
	1.5 Unit functions.....	15
	1.5.1 Cooling.....	15
	1.5.2 SuperFrost function.....	16
	1.5.3 Defrosting.....	16
	1.5.4 Temperature display.....	17
	1.5.5 Temperature memory.....	17
	1.5.6 Temperature monitoring with PC-KIT-STICK.....	17
	1.5.7 Temperature sensor.....	18
	1.6 Battery.....	18
	1.7 Additional unit components.....	18
	1.7.1 Temperature documentation.....	19
	1.7.2 Interior lighting.....	20
	1.7.3 Lockable glass door.....	20
2	Accessories	21
	2.1 GSM module.....	21
3	Safety	22
	3.1 Symbols in this instruction manual.....	22
	3.2 Purpose.....	23
	3.3 Foreseeable misuse.....	23
	3.4 Residual risks.....	24
	3.5 Safety markings.....	26
	3.6 Staff qualification.....	27
	3.7 Personal protective equipment.....	28
4	Transport and decommissioning	29
	4.1 Unit transport.....	29
	4.2 Final decommissioning of the unit.....	31
	4.3 Putting the unit back into operation.....	31
5	Set-up, installation and connection	32
	5.1 Setting up the unit.....	32
	5.2 Installing the unit.....	32
	5.3 Connecting the unit.....	33
6	Commissioning	34
	6.1 Activities during commissioning.....	34
	6.2 Programming the unit.....	34

6.2.1	Target temperature.....	35
6.2.1.1	Function of the target temperature.....	35
6.2.1.2	Displaying and changing the target temperature..	35
6.2.2	Temperature warning limits.....	37
6.2.2.1	Function of the temperature warning limits.....	37
6.2.2.2	Displaying and changing the temperature warning limits.....	38
6.2.3	Adjusting the humidity.....	40
7	Operation.....	41
7.1	Switching on the unit.....	41
7.2	Switching off the unit.....	41
7.3	Retrieving/erasing the temperature memory.....	42
7.4	Switching on the SuperFrost function.....	43
7.5	Setting up PC-KIT-NET (optional).....	43
7.6	Stocking the unit.....	48
8	Cleaning and disinfection.....	50
9	Maintenance.....	53
9.1	Safety inspection.....	53
10	Alarms.....	56
10.1	Alarm functions.....	56
10.2	Handling alarms.....	57
11	Status displays and error messages.....	58
11.1	Status displays.....	58
11.2	Error messages.....	58
12	Decommissioning and disposal.....	64
12.1	Decommissioning unit.....	64
12.2	Disposing of the unit.....	64
13	Appendix.....	66
13.1	Declaration of conformity.....	67
13.2	Technical data.....	69
13.3	Installation drawing.....	72
14	Index.....	74

1 Product description

1.1 Unit overview

Model LABO (example)



Fig. 1: Laboratory refrigerator (housing)

- 1 Key switch
- 2 Display and control unit ↪ *Chapter 1.2 “Display and control elements” on page 12*
- 3 Door lock
- 4 Door handle
- 5 Adjustable feet



Fig. 2: Laboratory refrigerator (interior)

- 1 Circulation cooling
- 2 Support rail
- 3 Supports
- 4 Drawer
- 5 Melt water container
- 6 Cooling machine

FROSTER LABO model (example)



Fig. 3: Laboratory freezer (housing)

- 1 Key switch
- 2 Display and control unit ↪ *Chapter 1.2 "Display and control elements" on page 12*
- 3 Door lock
- 4 Door handle
- 5 Adjustable feet



Fig. 4: Laboratory freezer (interior)

- 1 Circulation cooling
- 2 Support rail
- 3 Supports
- 4 Drawer
- 5 Melt water container
- 6 Cooling machine

1.2 Display and control elements

1.2.1 Design of the display and control unit

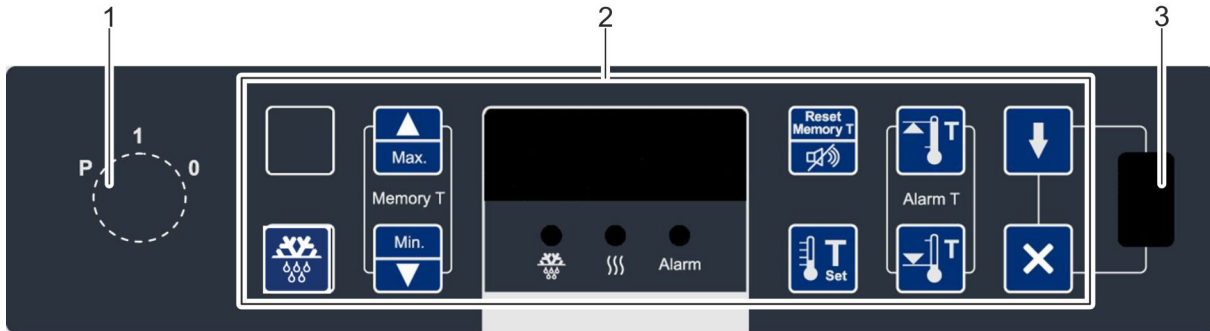


Fig. 5: Display and control unit

- 1 Key switch
- 2 Display and control unit (variable key assignment)
- 3 USB port

1.2.2 Function of the key switch

Tab. 1: Position of the key switch

Key position	Position	Function	Description
	"0"	Switch off unit.	<ul style="list-style-type: none"> ■ Switch off unit (for example to perform a restart). ■ The Display shows OFF. ■ To put the unit in a currentless state, unplug the power plug.
	"1"	Switch on unit.	<ul style="list-style-type: none"> ■ Operate unit. ■ Display and reset unit values (for example actual temperature or temperature limits). ■ Confirm alarms.
	"P"	Program unit.	<ul style="list-style-type: none"> ■ Adjust unit (for example target values or temperature warning limits). ■ Display parameter lists. ■ Confirm alarms.



Always operate unit in key position "1" in order to prevent manipulation. During normal operation, remove key and store it safely.

1.2.3 Function of buttons and displays

For the functions of buttons and displays, see the following overview:





Tab. 2: Buttons

Button	Description	Function
	[Max.]	Display maximum value of temperature memory.
	[Min.]	Display minimum value of temperature memory.
	[Reset]	Reset temperature memory. Switch off buzzer.
	[Temperature warning max.]	Read upper temperature warning limit.
	[Temperature warning min.]	Read lower temperature warning limit.
	[Target temperature]	Read target temperature.
	[Defrosting]	Activate additional defrosting.
	[Additional light / additional defrosting]	Combination button: Activate additional defrosting. Switch interior lighting on/off permanently.
	[SuperFrost]	Activate SuperFrost function.
	[Humidity]	Adjust humidity.
	[Light]	Switch interior lighting on/off permanently.
	[Start]	Start data read-out via a PC-KIT-STICK.
	[Cancel]	Cancel data read-out via a PC-KIT-STICK.
	[Unassigned button]	Unassigned button without function.



Depending on the version, certain buttons have a multiple function in combination with other buttons.

Tab. 3: Displays

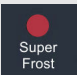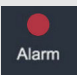Display	Description	Function
	"Defrosting"	Defrosting is active.
	"Humidity"	Humidity is active (temperature consistency improved, humidity high).
	"SuperFrost"	SuperFrost is active.
	"Alarm"	An alarm was triggered.

1.3 Scope of delivery

Interior fittings

The following interior fittings approved for operation are supplied for each unit in line with the unit specifications:

- Wire shelves
- Drawers

Lockable unit door

The unit is fitted with a lockable unit door.

Keys included in delivery

Depending on the unit components, the following keys are included in delivery:

- Unit key
- Door key (for lockable unit door)

1.4 Interfaces

The unit is equipped with the following interface for connecting additional units (hereafter referred to as "modules") for monitoring and documenting the temperature:

Tab. 4: Interfaces

Interface	Module
LAN interface	PC-KIT-NET (optional)
USB port	PC-KIT-STICK
Potential-free alarm contact	Remote warning system (for example GSM-MODUL or connection to building control system (see circuit diagram on unit))

For information about connecting the modules to the unit, see the relevant product documentation.



Connect only compatible USB flash drives!

Connect only modules with the following properties to the USB port:

- Maximum memory capacity of 32 GB
- Formatted in FAT-32 format



Remove USB flash drive after use!

The USB flash drive must not remain permanently in the unit.

- Remove USB flash drive after the data transfer is completed → “Error messages and status displays of the PC-KIT-STICK” on page 63.

1.5 Unit functions

1.5.1 Cooling

Cooling machine

The unit is equipped with a cooling machine with interior evaporator for cooling chilled goods.

After initial commissioning or recommissioning, the cooling machine takes a while to cool the interior down to the set target temperature value.

Circulation cooling



With circulation cooling, the air circulates around the interior of the refrigerator. Circulation cooling is performed by means of ventilation plates and a fan.

Circulation cooling reduces physically induced temperature differences and the target temperature is kept constant throughout the interior.

Circulation cooling is set to continuous operation and is switched off automatically when the door is opened.

Condenser



The condenser transfers the generated thermal energy to the ambient air.



Depending on the model, the condenser is fitted either on the back wall or in the machine room.

Fig. 6: Condenser (example)

1.5.2 SuperFrost function



The unit is equipped with a SuperFrost function.

With the SuperFrost function, it is possible to cool the unit down to the lowest possible temperature in the shortest time ↪ [Chapter 7.4 “Switching on the SuperFrost function” on page 43.](#)

1.5.3 Defrosting

Automatic defrosting



The unit defrosts automatically every 12 hours.

Automatic defrosting is time and temperature-controlled.

During automatic defrosting, the system ensures that the unit maintains the target temperature value.

Melt water is conducted into the exterior melt water container and evaporates there.



When the defrosting process is active, the “Defrosting” LED display lights up on the Display.

Additional defrosting

In addition to automatic defrosting, it is possible to start the defrosting process manually.



➔ Press button [Additional defrosting] for four seconds.

⇒ The defrosting process begins.

1.5.4 Temperature display



Fig. 7: Display

The temperature is shown on the Display of the unit's display and control unit ↪ *Chapter 1.2 "Display and control elements" on page 12.*

The Display indicates the temperature of the chilled goods.

i *The temperature display does not indicate the current air temperature of the interior.*

The alarm function and the temperature warning limits are based on the temperature of the chilled goods as shown on the Display.

1.5.5 Temperature memory




The temperature memory records the maximum and minimum temperature values reached during operation.

The temperature memory begins recording when the target temperature is reached or, at the latest, after two hours of operation.

The temperature memory is reset automatically when the unit is restarted or reset.

The temperature memory can be retrieved and erased manually ↪ *Chapter 7.3 "Retrieving/erasing the temperature memory" on page 42.*

1.5.6 Temperature monitoring with PC-KIT-STICK

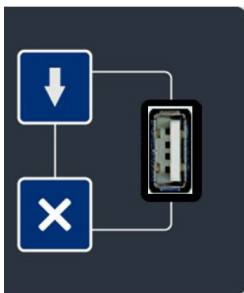

The temperature memory of the unit documents the maximum and minimum temperature values.

PC-KIT-STICK is the easiest method of automatic temperature documentation.




The data transfer is performed via a USB flash drive. Even while the data are read out, the temperatures is constantly documented. An unlimited number of units can be incorporated ↪ *"Interfaces" on page 15.*

For constant temperature monitoring, the unit can be equipped with temperature documentation via the network ↪ *Chapter 1.7.1 "Temperature documentation" on page 19.*

1.5.7 Temperature sensor

The unit is equipped with multiple temperature sensors.

For the functions of the respective temperature sensors, see the following table:

Colour	Description	Function
	Evaporator sensor (red)	Controls the defrosting phase
	Control sensor (brown)	Controls the monitoring unit
	Control sensor (violet)	Controls the monitoring unit

1.6 Battery

Battery function

The battery powers the temperature monitoring for up to 30 hours in the event of a power failure.

Temperature monitoring ensures that the temperature progression is stored and the temperature alarm is triggered if necessary.

The battery charges during normal operation and is monitored by an automatic charging system.

Technical data for the battery:

- 6 V, 4 Ah



The battery does not supply the unit!

The battery only provides power to the temperature monitoring. The battery does not provide a back-up to the unit and therefore ensure that the interior temperature is maintained.

1.7 Additional unit components

The unit is equipped with standard equipment.

The standard equipment can be supplemented with additional unit components.

The following additional components are available for the unit:


- Temperature documentation via the network
- Interior lighting

1.7.1 Temperature documentation





The unit can optionally be equipped with additional temperature documentation.

The following table contains an overview of additional unit components.

For more information on installation and use, see the documentation for the relevant product.

For information about the software for temperature documentation, see the  "Assembly and connection manual – PC-KIT-NET".

Tab. 5: Temperature documentation options

Description	Function	Figure
PC-KIT-NET	<p>Automatic temperature documentation and monitoring via the network (LAN).</p> <p>Unlimited connection of units possible.</p> <p>Simultaneous data access for up to 20 clients per server. Two server accesses possible.</p>	
Disc-type pen-recording thermometer	<p>Depending on the model, it is either installed in the machine room cover of the refrigerator or freezer or integrated in an additional housing extension.</p>	
Pen-recording thermometer for placing in the unit	<p>The easiest way to retrofit temperature documentation.</p> <p>Temperature measuring range from +25 °C to -40 °C.</p> <p>Not suitable for FROSTER BL 180 PRO-ACTIVE and FROSTER BL 650 PRO-ACTIVE</p>	
External temperature documentation	<p>The sensor makes it possible to measure the temperature and transfer the data to existing interfaces on site.</p> <p>The user/owner may have to check for compatibility.</p>	

1.7.2 Interior lighting

Interior lighting (optional)

The unit is optionally equipped with interior lighting.

The interior lighting automatically switches on when the door is opened and switches off when the door is closed.



In units with a glass door, the interior lighting can be permanently switched on or off with the *[Light]* button.

1.7.3 Lockable glass door

Door lock

Your unit can be equipped with a lockable glass door.

For more information on installing the glass door, contact KIRSCH.

2 Accessories

The unit can be equipped with the following accessories:

- GSM module

2.1 GSM module



Fig. 8: GSM module

The GSM module is used to forward alarm messages to a mobile or fixed-line telephone network via text message.



To operate the GSM module you need a SIM card (not supplied).

A GSM module can manage up to three units. You can use the supplied software to adjust the alarm limits (upper and lower warning limit). You need a temperature sensor to do so (not supplied).

The GSM module has three inputs:

- Two inputs for the potential-free contact
- One input for the optional temperature sensor

Alarm messages are forwarded to the GSM module via the potential-free contact (connection for remote warning system). When the alarm limits are reached, an optical and acoustic signal is triggered on the GSM module and an alarm text message is sent.

The acoustic alarm is acknowledged by sending a confirmation text message to the GSM module or by pressing a button on the GSM module. The optical alarm remains until the fault has been rectified and also acknowledged.

Units without a potential-free contact can be retrofitted with a temperature sensor.

3 Safety

This section provides an overview of all important safety aspects for optimal protection of patients and staff, and for safe and trouble-free operation of the unit.

Non-compliance with the instructions and safety warnings in this instruction manual can cause considerable risks.

3.1 Symbols in this instruction manual

Safety instructions

Safety instructions are marked with symbols in this instruction manual. Safety instructions are initiated by signal words that express the degree of risk.

In order to avoid accidents, injury and damage and ensure maximum patient safety, always comply with safety instructions and act with care.



DANGER!

This combination of symbol and signal word indicates an immediately hazardous situation that will lead to death or serious injury unless avoided.



WARNING!

This combination of symbol and signal word indicates a potentially hazardous situation that can lead to death or serious injury unless avoided.



CAUTION!

This combination of symbol and signal word indicates a potentially hazardous situation that can lead to minor or slight injury unless avoided.



NOTICE!






This combination of symbol and signal word indicates a potentially hazardous situation that can lead to property damage or environmental damage unless avoided.

Hints and recommendations

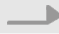

This symbol highlights useful hints and recommendations as well as information for efficient and trouble-free use of the unit.

Other markings

Mark	Explanation
	Step-by-step instructions
	Results of actions
	References to sections in this instruction manual
	Lists without a specified order
	References to the instruction manuals for accessories and optional parts

3.2 Purpose

The laboratory refrigerator and laboratory freezer are used for commercial storage of biological and chemical substances. The laboratory refrigerator and the laboratory freezer are not medical products.

3.3 Foreseeable misuse

The unit is not designed for domestic use. The unit is used for commercial storage of chilled goods in line with its intended purpose.

Do not use the unit to cool warm goods. Do not store chilled goods in the unit if their cooling chain was interrupted during delivery or stock transfer.

Do not store food or drink in the unit.

Do not store chilled goods that exceed the carrying capacity of the wire shelves and drawers.

3.4 Residual risks

Infection of the user



WARNING!

Danger of infection due to inadequate hygiene, disinfection and sterilisation!

Contact with parts that have not been cleaned, sterilised or disinfected poses the danger of infection.

- Comply with the applicable requirements for personal hygiene, disinfection and sterilisation.
- Clean, disinfect and sterilise the unit before filling it for the first time.
- Clean, disinfect and sterilise the unit in line with this instruction manual during operation.

Escaping refrigerant



WARNING!

Danger due to escaping refrigerant!

The refrigerant used in the unit is explosive and pressurised, and can cause serious injury if it comes into contact with the eyes or skin.

- During transport and set-up, do not bend or pierce the tubing and the evaporator.
- Do not damage the surface coat (scratching it off, for example).
- Wear safety glasses and protective gloves when handling the refrigerant circuit.

Hot surface



WARNING!

Danger due to hot surface!

The marked areas of the unit can cause severe skin injuries if touched.

- Do not touch areas of the unit marked in this way.
- These areas are very hot and can still cause burns several hours after the unit has been switched off.

Eco cooling machines



WARNING!

Danger due to unauthorised intervention in the cooling circuit!

The cooling machine contains the natural refrigerant propane R290 / isobutane R600a. The refrigerant used in the unit is flammable and can cause serious injury if it comes into contact with the eyes or skin. Unauthorised intervention in the cooling circuit poses a danger of injury!

- On site, only the entire cooling machine can be replaced.

The cooling machine is a permanently technically sealed system in line with EN 1127-1.

Contamination of chilled goods



CAUTION!

Danger to chilled goods due to contamination!

A dirty unit can contaminate the stored chilled goods. A contamination can damage or destroy the chilled goods.

- Do not store contaminated chilled goods in the unit.
- After any contamination, clean, disinfect and sterilise the unit and the chilled goods.
- Clean, disinfect and sterilise the unit in line with this instruction manual during operation.
- Comply with the applicable requirements for personal hygiene, disinfection and sterilisation.

Blocked interior ventilation



CAUTION!

Danger of damage and injury due to blocked interior ventilation!

Without adequate ventilation of the cooling machine the air circulation inside the unit is no longer ensured. This can lead to a temperature drop inside the unit, which can damage the chilled goods.

- Make sure that the ventilation in the upper area is clear.
- Do not cover ventilation grille with chilled goods.

Missing power supply



NOTICE!

When a power failure occurs, the cooling unit is switched off, which can damage or destroy the chilled goods.

- Make sure that the unit is protected from power failures (for example by an uninterruptible power supply).
- After a power failure, make sure that the cooling unit is supplied with power again.
- If necessary move the chilled goods to another location and inform the service department.

Blocked outside ventilation



NOTICE!

Danger of damage due to blocked outside ventilation!

Without adequate ventilation of the cooling machine the unit can overheat and be damaged.

- Do not set up units in close proximity to each other.
- Make sure that the ventilation (back and front) on the unit is clear.
- Do not cover the cooling machine.
- Allow only qualified staff to install the unit.

3.5 Safety markings

Type plate

Modell Type	BL 176 PRO-ACTIVE	Kühl-Inhalt Gross volume	170 Liter
Fabrik-Nr. Serial-No.	170 30 25000	Netto-Inhalt Net volume	150 Liter
Kältemittel Refrigerant	R600a 0,030 kg	Kühl-Aggregat Cooling system	VKS 2716 KS
Wechselstrom Voltage	220-240V~ 50 Hz	Stromaufnahme Power input	0,53 A
Klimaklasse FD	0,1,2,3,4,5,6,7	Zusatzheizungen mit Leistungen > 100W	---

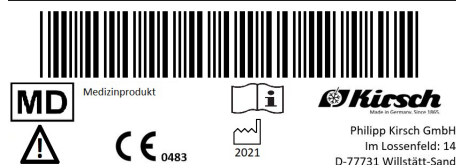


Fig. 9: Example type plate (BL 176 PRO-ACTIVE)

The type plate is mounted in the interior of the unit.

The type plate contains the following information:

- Model
- Factory no.
- Refrigerant type and quantity
- Alternating current
- Climate class
- Refrigerator capacity (volume of entire interior)
- Net capacity (volume of usable interior)
- Cooling unit
- Power consumption
- Additional heaters with output > 100 W

3.6 Staff qualification

Inadequate staff qualification



WARNING!

Danger of damage and injury due to commissioning by unqualified staff!

If commissioning is performed incorrectly by unqualified staff, serious damage to the chilled goods can result, which in turn can lead to serious injuries.

- Have all tasks performed only by staff qualified for those tasks.
- Keep unauthorised persons away from the working area.

Staff qualifications

This manual specifies the staff qualifications for various fields of activity as listed below:

System/network administrator (recommended)

The system/network administrator has the training, IT skills and experience required to set up the system requirements and working environment, including all technical equipment, to enable the software to be used.

The system/network administrator performs the following duties:

- Installing KIRSCH-DATANET
- Integrating the unit in the network

The system/network administrator has been authorised by the owner to manage the users of the software and to make settings to the software.

Unit officer

The unit officer is the person who is nominated for this task by the owner of the unit and has received instruction on their duties.

Enter the name of the unit officer and the date of instruction in the medical product book and confirm with your signature.

The unit officer meets the following requirements:

- The unit officer knows the intended purpose, the foreseeable misuse and the residual risks of the unit.
- The unit officer is familiar with the instruction manual and all other safety-related documents.
- The unit officer has been instructed in the technically correct and safe handling of the unit.

The unit officer performs the following tasks:

- The unit officer instructs users in how to handle the unit.

User

The user is the person who uses and operates the unit according to its intended purpose. The unit may only be used and operated by trained specialist staff.

The user has been instructed in the technically correct and safe use of the unit in accordance with the relevant laws and ordinances.

General staff qualification requirements

Staff members must be persons who can be expected to perform their work reliably. Persons whose reactions are impaired, e.g. by drugs, alcohol or medication, are not permitted.

When choosing employees, observe the age and vocation regulations that apply at the deployment site.

3.7 Personal protective equipment

Personal protective equipment protects staff members from dangers that could affect their safety or health at work.

Always wear the personal protective equipment specified in the various chapters of this manual before starting the relevant task.

Observe instructions on personal protective equipment that are installed in the working area.

When performing various tasks on and with the unit, staff must wear personal protective equipment. This is indicated specifically in the individual chapters in this manual. This personal protective equipment is explained below:



Chemical-resistant protective gloves

Chemical-resistant protective gloves protect the hands from aggressive chemicals.



Protective gloves

Protective gloves protect the hands from friction, abrasions, piercing or deeper injuries as well as from contact with hot surfaces.



Safety shoes

Safety shoes protect the feet from crushing, from falling parts, and from slipping on slippery ground.

4 Transport and decommissioning

4.1 Unit transport

On delivery, the unit is transported after consultation with the specialist dealer.

When transporting the unit during a change of location while the unit is still operating, observe the following safety instructions.

Safety during transport



WARNING!

Danger of crushing injuries from falling unit!

If the unit is tilted, it will fall over in an uncontrolled fashion. When the unit falls over, there is a danger of crushing to the hands and feet.

- Wear safety shoes and safety gloves when transporting the unit.
- Transport the unit in an upright position.



WARNING!

Danger due to escaping refrigerant!

The refrigerant used in the unit is explosive and pressurised, and can cause serious injury if it comes into contact with the eyes or skin.

- During transport and set-up, do not bend or pierce the tubing and the evaporator.
- Do not damage the surface coat (scratching it off, for example).
- Wear safety glasses and protective gloves when handling the refrigerant circuit.



NOTICE!

Danger of damage to the melt water container!

The units in the table below are equipped with a melt water container on their underside, which can be damaged by incorrect transport:

- Place unit on pallet and transport it.
- Do not lift unit with forklift or pallet truck unless it is on a pallet.

LABO 520 PRO-ACTIVE	LABO 520 ULTIMATE
LABO 720 PRO-ACTIVE	LABO 720 ULTIMATE
LABO 520 CR PRO-ACTIVE	LABO 520 ULTIMATE
LABO 720 CR PRO-ACTIVE	LABO 720 ULTIMATE
LABO 720 CHROMAT PRO-ACTIVE	LABO 720 CHROMAT ULTIMATE
FROSTER LABO 530 PRO-ACTIVE	FROSTER LABO 530 ULTIMATE
FROSTER LABO 730 PRO-ACTIVE	FROSTER LABO 730 ULTIMATE
FROSTER LABO 530 CR PRO-ACTIVE	FROSTER LABO 530 CR ULTIMATE
FROSTER LABO 730 CR PRO-ACTIVE	FROSTER LABO 730 CR ULTIMATE

Transporting the unit to a new location

- Protective equipment: ■ Safety shoes
 ■ Protective gloves

Decommission the unit as shown in section ↗ “Final decommissioning of the unit” on page 31.

- 1.** ➡ Transport the unit to the new location.




Waiting time before recommissioning:

Once the unit has been transported, set up the unit and wait one hour before recommissioning, so that the refrigerant can distribute itself evenly throughout the refrigerant circuit again.

- 2.** ➡ Recommission the unit as shown in section ↗ “Putting the unit back into operation” on page 31.

4.2 Final decommissioning of the unit


Final decommissioning of the unit Insert unit key in key switch.

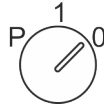

1. ➤ Turn unit key to position "0".
2. ➤ Unplug power plug.
3. ➤ Remove unit key.
4. ➤ Leave unit door open to prevent the formation of odours and mould.

4.3 Putting the unit back into operation

Putting the unit back into operation Before recommissioning, clean and disinfect the unit ↪ *Chapter 8 "Cleaning and disinfection" on page 50.*

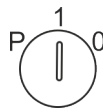

1. ➤ Plug in power plug.
2. ➤ Insert unit key in key switch.
3. ➤ Turn unit key to position "1".
⇒ The Display shows the interior temperature.
4. ➤ Check target temperature and adjust if necessary ↪ *Chapter 6.2.1.2 "Displaying and changing the target temperature" on page 35.*



Lead time of the temperature alarm!

The temperature alarm is active at the earliest when the target temperature is reached, and at the latest after two hours of unit operation.

5. ➤ When the unit has reached its target temperature, place the chilled goods in the unit.

5 Set-up, installation and connection

5.1 Setting up the unit

Set-up

When setting up the unit, observe the following set-up conditions:

- Ensure an ambient temperature between +10 °C – +38 °C (as per climate class SN and ST).
- Ensure an ambient temperature between +10 °C – +32 °C (as per climate class SN).
- Select a dry and well ventilated set-up location.
- Do not set up the unit next to heat sources.
- Avoid direct sunlight.
- Set up the unit on a firm, level surface.
- Compensate for uneven ground with adjustable feet.
- Do not cover, block or line the ventilation grille.
- Do not cover, block or line the cooling machine.
- Maintain a minimum distance of 2 cm between wall and unit (for example by means of a spacer).
- Maintain a minimum distance of 2 cm between units.
- Comply with the minimum available space at the set-up location in order to avoid a potentially explosive as defined in DIN EN378-1 for the use of refrigerators and freezers with combustible refrigerants (R600a, R290a).



The basis for calculating the minimum available space at the set-up location is the largest fill quantity of refrigerant in the existing or planned refrigerators. For the respective fill quantity of refrigerant, see table ↗ Chapter 13.2 “Technical data” on page 69.

R600a: volume = refrigerant fill quantity / 8.6 g/m³

R290a: volume = refrigerant fill quantity / 7.6 g/m³

5.2 Installing the unit

Overview of units suitable for installation

LABO 100 PRO-ACTIVE

LABO 126 PRO-ACTIVE

Installing

Have the unit installed only by staff qualified for the task. If in doubt, contact KIRSCH ↗ “Manufacturer’s address” on page 4.

Install the unit as shown in the installation drawing provided.

5.3 Connecting the unit

Connecting



The unit is designed according to protection class I and protection type IP 20 and is ready to be plugged in.

Ensure the connection conditions specified in the technical data and on the type plate.

1. ➤ Check the connecting cable of the power plug for damage.
2. ➤ Plug in power plug.

6 Commissioning

Personnel: ■ Unit officer

6.1 Activities during commissioning

Commissioning sequence



WARNING!

Danger of damage and injury due to commissioning by unqualified staff!

If commissioning is performed incorrectly by unqualified staff, serious damage to the chilled goods can result, which in turn can lead to serious injuries.

- Have all tasks performed only by staff qualified for those tasks.
- Keep unauthorised persons away from the working area.

Commissioning consists of the following activities:

- 1 - Clean and disinfect the interior of the unit ↪ *Chapter 8 “Cleaning and disinfection” on page 50.*
- 2 - Switch on the unit ↪ *Chapter 7.1 “Switching on the unit” on page 41.*
- 3 - Program the unit ↪ *Chapter 6.2 “Programming the unit” on page 34.*
- 4 - Wait until the target temperature is reached.
- 5 - Stock the unit ↪ *Chapter 7.6 “Stocking the unit” on page 48.*

6.2 Programming the unit

Ensure only authorised personnel have access to the key



NOTICE!

Danger due to faulty programming!

By putting the key in key position “P” it is possible to access the programming functions of the unit. Programming that is unsuitable for the chilled goods can cause damage to the chilled goods.

- Have programming performed by qualified employees.
- Once programming is complete, turn the key switch to position “1” and remove the key.
- Secure the key against unauthorised access.
- Do not operate the unit in key position “P”.

6.2.1 Target temperature

6.2.1.1 Function of the target temperature

The target temperature specifies the temperature at which the unit is operated to store the chilled goods in optimal conditions.

The target temperature of the unit is preset by KIRSCH.



Changes to the target value do not change the temperature warning limits. These are adjusted manually ↪ Chapter 6.2.1.2 “Displaying and changing the target temperature” on page 35.

6.2.1.2 Displaying and changing the target temperature

Changes must be performed by qualified employees



NOTICE!

The factory-set values for the target temperature comply with the valid DIN standards. Incorrect changes to the target temperature and the temperature warning limits can cause irreparable damage to the chilled goods.

Displaying the target temperature



→ Press button [Target temperature].

⇒ The current target temperature is displayed.

Changing the target temperature

1. ➤ Insert unit key in key switch.




NOTICE!

Danger due to faulty programming!

By putting the key in key position “P” it is possible to access the programming functions of the unit. Programming that is unsuitable for the chilled goods can cause damage to the chilled goods.

- Have programming performed by qualified employees.
- Once programming is complete, turn the key switch to position “1” and remove the key.
- Secure the key against unauthorised access.
- Do not operate the unit in key position “P”.

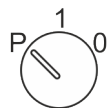


2. ➤ Turn unit key to position “P”.



3. ➤ Press and hold button [*Target temperature*].



4. ➤ Use button [*Min.*] or [*Max.*] to set the required target value.

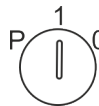

5. ➤ Turn unit key to position “1”.

⇒ The unit regulates the temperature to the specified target value.



Monitor the temperature progression!

The unit will not reach the new target value immediately after the change is made.

Monitor the temperature progression on the Display or using optional temperature documentation (e.g. KIRSCH-PC-KIT).

6. ➤ Read the current temperature and monitor the subsequent temperature progression until the target temperature is reached.



Consequences of modified target temperature

As soon as the target temperature is reached, the temperature warning limits must be adjusted so that the target temperature is above or below the temperature warning limits. Otherwise the temperature alarm will be triggered ↪ Chapter 6.2.2.2 “Displaying and changing the temperature warning limits” on page 38.

The set target temperature is saved automatically. The target temperature is retained after a power failure and after the unit has been switched off.

7. ➔ When the unit has reached its target temperature, stock the unit.

6.2.2 Temperature warning limits

6.2.2.1 Function of the temperature warning limits

The temperature warning limits define how much deviation the unit will tolerate between the actual temperature and the target temperature. The temperature warning limits are unit-specific. They can be adapted to the requirements of the chilled goods.

The values for the temperature warning limits are set at the factory and comply with the valid DIN standards for the unit.



Recommended temperature warning limits

The temperature warning limits must not be the same as the target temperature.

Set the temperature limits as follows:

- *Upper temperature warning limit: at least 3 °C higher than the target temperature*
- *Lower temperature warning limit: at least 3 °C lower than the target temperature*



Recommended temperature warning limits for FROSTER

To ensure that the unit functions properly, the temperature warning limits must be set for FROSTER models as shown in table ↪ “Temperature warning limits” on page 38.

Temperature warning limits

Tab. 6: Temperature warning limits for LABO and FROSTER LABO

Model	Lower temperature warning limit	Target temperature	Upper temperature warning limit
LABO 100 PRO-ACTIVE LABO 288 PRO-ACTIVE LABO 340 PRO-ACTIVE LABO 468 PRO-ACTIVE LABO 520 PRO-ACTIVE LABO 720 PRO-ACTIVE LABO 720 CHROMAT PRO-ACTIVE LABO 288 ULTIMATE LABO 340 ULTIMATE LABO 468 ULTIMATE LABO 520 ULTIMATE LABO 720 ULTIMATE LABO 720 CHROMAT ULTIMATE	+2 °C	+5 °C	+8 °C
FROSTER LABO 330 PRO-ACTIVE FROSTER LABO 530 PRO-ACTIVE FROSTER LABO 730 PRO-ACTIVE FROSTER LABO 330 ULTIMATE FROSTER LABO 530 ULTIMATE FROSTER LABO 730 ULTIMATE	-55 °C	-25 °C	-15 °C

6.2.2.2 Displaying and changing the temperature warning limits

Displaying the temperature warning limits

Upper temperature warning limit



1. → Press button [Temperature warning max.].
⇒ The upper temperature warning limit is shown on the Display.

Lower temperature warning limit



2. → Press button [Temperature warning min.].
⇒ The lower temperature warning limit is shown on the Display.

Changing the temperature warning limits



Recommended temperature warning limits

The temperature warning limits must not be the same as the target temperature.

Set the temperature limits as follows:

- Upper temperature warning limit: at least 3 °C higher than the target temperature
- Lower temperature warning limit: at least 3 °C lower than the target temperature




Recommended temperature warning limits for FROSTER

To ensure that the unit functions properly, the temperature warning limits must be set for FROSTER models as shown in table ↗ “Temperature warning limits” on page 38.

1. ➔ Insert unit key in key switch.

2. ➔ Turn unit key to position “P”.

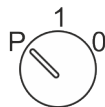

3. ➔ Press and hold button [Temperature warning max.].




4. ➔ Press button [Max.] repeatedly until the value for the new upper temperature warning limit is reached.



5. ➔ Press button [Min.] repeatedly until the value for the new lower temperature warning limit is reached.



6. ➔ Turn unit key to position “1”.

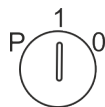

⇒ The Display shows the current interior temperature.



The set temperature warning limits are saved automatically. The temperature warning limits are retained after a power failure and after the unit has been switched off.

6.2.3 Adjusting the humidity

Overview of units with adjustable humidity

LABO 288 PRO-ACTIVE	LABO 288 ULTIMATE
LABO 340 PRO-ACTIVE	LABO 340 ULTIMATE
LABO 468 PRO-ACTIVE	LABO 468 ULTIMATE
LABO 520 PRO-ACTIVE	LABO 520 ULTIMATE
LABO 520 CR PRO-ACTIVE	LABO 520 ULTIMATE
LABO 720 PRO-ACTIVE	LABO 720 ULTIMATE
LABO 720-CHROMAT PRO-ACTIVE	LABO 720-CHROMAT ULTIMATE

Adjusting the humidity

Two humidity levels are available for the interior of the unit:


high humidity - approx. 90%

low humidity - approx. 65%



If the “high humidity” setting is selected, the recirculating blower in the interior is constantly in operation and is switched off only when the door is opened. If the “low humidity” setting is selected, the blower switches off briefly each time the cooling machine stops.

1. ▶ Insert the unit key in the key switch.

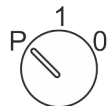


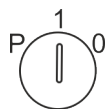
2. ▶ Turn unit key to position “P”.



3. ▶ Press button [Humidity] and hold it for at least four bseconds.

4. ▶ Select humidity.

⇒ After you select “high humidity”, LED display  lights up.
After you select “low humidity”, the LED display goes out.




5. ▶ Turn unit key to position “1”.

7 Operation

Personnel: ■ User

7.1 Switching on the unit

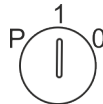

1. ➤ Insert the key into the key switch.
2. ➤ Turn unit key to position "1".
⇒ The Display shows the interior temperature.
3. ➤ Remove the key and store it so it is protected against unauthorised access.
4. ➤ Check target temperature and adjust if necessary ↪ *Chapter 6.2.1.2 "Displaying and changing the target temperature" on page 35.*



Lead time of the temperature alarm!

The temperature alarm is active at the earliest when the target temperature is reached, and at the latest after two hours of unit operation.


5. ➤ Stock the unit when the target temperature is reached.



The door is stiff when it is first opened

Cooling creates a vacuum in the interior, so you may require additional force when opening the door for the first time.

7.2 Switching off the unit

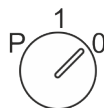


1. ➤ Insert unit key in key switch.
2. ➤ Turn unit key to position "0".
⇒ The unit switches to standby mode and the Display shows 0FF.



Only switch the stocked unit off for short periods!

To protect the chilled goods, only switch the stocked unit off for short periods.

To shut down the unit for a longer period, proceed as described in  Chapter 4.2 "Final decommissioning of the unit" on page 31.

7.3 Retrieving/erasing the temperature memory

Retrieving the temperature memory



The temperature memory records the temperature from the moment the target temperature is reached.

Restarting the unit resets the temperature memory.



1. Press button [Max.].

⇒ The maximum value of the temperature memory is displayed.



2. Press button [Min.].

⇒ The minimum value of the temperature memory is displayed.

Erasing the temperature memory



1. Press and hold button [Max.].



2. Press button [Reset].

⇒ The Display shows .

The maximum value of the temperature memory is erased.



3. Press and hold button [Min.].



4. Press button [Reset].

⇒ The Display shows .


The minimum value of the temperature memory is erased.

7.4 Switching on the SuperFrost function

The lower temperature warning limit of the unit is set to the lowest value ↪ Chapter 6.2.2 “Temperature warning limits” on page 37.

Switching on SuperFrost



1. ➤ Press button [SuperFrost].
⇒ LED display  lights up.



The SuperFrost function runs for 12 hours and ends automatically.

Manually ending SuperFrost



2. ➤ Press button [SuperFrost] again.
⇒ The SuperFrost function is deactivated.
3. ➤ Adapt lower temperature warning limit to the target temperature ↪ Chapter 6.2.2 “Temperature warning limits” on page 37.

7.5 Setting up PC-KIT-NET (optional)

Tab. 7: Factory settings TCP/IP

IP address of unit	192.168.0.101
Subnet mask	255.255.255.0
IP address of standard gateway	192.168.0.200



If you want to commission more than one unit, the following applies:

Always connect one unit after the other to the network and complete the network connection on the PC.



The IP address of the unit and the IP address of the standard gateway must be configured in the same address field:

IP address of unit (example)	<u>180.160.15.1</u>
IP address of standard gateway (example)	<u>180.160.15.2</u>





For temperature monitoring via the network, the KIRSCH-DATANET software (version 5.0 or higher) must be installed (complete assembly and connection manual is supplied with PC-KIT-NET).

Setting up the IP address on the unit



Every unit needs its own IP address that is not yet in use in your local network.

At the factory, the IP address is set to:
192.168.0.101

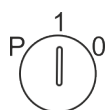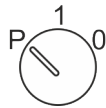


1. ➤ Insert unit key in key switch.
2. ➤ Set unit key to "P".
3. ➤ Press and hold [Max.] and [Min.] simultaneously for four seconds.
⇒ The display shows .
4. ➤ Press [Max.] or [Min.] to go to level and confirm with Set.
5. ➤ Press [Min.] to select the parameter .
6. ➤ Press and hold [Target temperature].
⇒ The first three digits of the IP address are displayed.
7. ➤ Use [Max.] and [Min.] to set the selected digits of the required IP address.
⇒ The number changes on the display in steps of one.
8. ➤ Release [Target temperature].
⇒ The first three digits of the IP address have been entered.
9. ➤ Press [Min.] to select the next parameter .
10. ➤ Press and hold [Target temperature].
⇒ The next digits of the IP address are displayed.
11. ➤ At levels L52 and L53 repeat steps 6 to 10 until the IP address of the unit has been entered completely.
12. ➤ Set key switch to "1" to complete your input.
⇒ The IP address of the unit has been set up.

Setting up the subnet mask on the unit



At the factory, the subnet mask is: 255.255.255.0

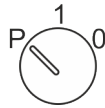

1. Insert unit key in key switch.

2. Set unit key to "P".



3. Press and hold [Max.] and [Min.] simultaneously for four seconds.

⇒ The display shows **Rd.r.**



4. Press [Max.] and [Min.] to go to level **L5r** and confirm with Set.



5. Press [Min.] to select the parameter **L55**.



6. Press and hold [Target temperature].

⇒ The first three digits of the subnet mask are displayed.



7. Use [Max.] and [Min.] to set the selected digits of the required subnet mask.

⇒ The number changes on the display in steps of one.



8. Release [Target temperature].

⇒ The first three digits of the subnet mask have been entered.




9. Press [Min.] to select the next parameter **L56**.

10. Press and hold [Target temperature].

⇒ The next digits of the subnet mask are displayed.

11. At levels L57 and L58 repeat steps 6 to 10 until the subnet mask has been entered completely.

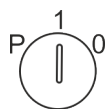

12. Set key switch to "1" to complete your input.



⇒ The subnet mask has been set up.

Setting up the IP address of the standard gateway on the unit



The standard gateway of each unit needs its own IP address that is not yet in use in your local network.

At the factory, the IP address of the standard gateway is set to: 192.168.0.200

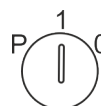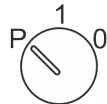

1. ➤ Insert unit key in key switch.
2. ➤ Set unit key to "P".
3. ➤ Press and hold [Max.] and [Min.] simultaneously for four seconds.
⇒ The display shows .
4. ➤ Press [Max.] and [Min.] to go to level and confirm with Set.
5. ➤ Press [Min.] to select the parameter .
6. ➤ Press and hold [Target temperature].
⇒ The first three digits of the IP address are displayed.
7. ➤ Use [Max.] and [Min.] to set the selected digits of the required IP address.
⇒ The number changes on the display in steps of one.
8. ➤ Release [Target temperature].
⇒ The first three digits of the IP address have been entered.
9. ➤ Press [Min.] to select the next parameter .
10. ➤ Press and hold [Target temperature].
⇒ The next digits of the IP address are displayed.
11. ➤ At levels L62 and L63 repeat steps 6 to 10 until the IP address of the standard gateway has been entered completely.
12. ➤ Set key switch to "1" to complete your input.
⇒ The IP address of the standard gateway has been set up.

Connecting the unit to the local network


- Personnel: ■ System/network administrator (recommended)
- Materials: ■ Network cable (EIA/TIA-568 standard), (included in the scope of delivery of PC-KIT-NET)

Requirements:

- The TCP/IP module has been installed (↪ *Chapter 1.7 “Additional unit components” on page 18* and ↪ *“Assembly and connection manual – PC-KIT-STICK/PC-KIT-NET”*).
- The KIRSCH-DATANET software (version 5.0 or higher) has been installed on the local PC or local network.

1. ➤ Switch off the unit ↪ *Chapter 7.2 “Switching off the unit” on page 41*.
2. ➤ Unplug the power plug.
3. ➤ Remove the dust cover from TCP/IP output.
4. ➤ Insert the network cable in the TCP/IP output.
5. ➤ Connect the network cable to the network socket.
⇒ The unit has been connected to the network.
6. ➤ Insert the power plug in the socket.
7. ➤ Switch on the unit ↪ *Chapter 7.1 “Switching on the unit” on page 41*.

Reading out the MAC IP address on the unit

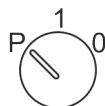

1. ➤ Insert unit key in key switch.
2. ➤ Set unit key to “P”.
3. ➤ Press and hold *[Max.]* and *[Min.]* simultaneously for four seconds.
⇒ The display shows .
4. ➤ Press *[Max.]* and *[Min.]* to go to level and confirm with *Set*.
5. ➤ Press *[Min.]* to select the parameter .
- ⇒ The first digits of the MAC address are displayed.




The first two displayed digits “00” are not part of the MAC address.



6. Press [Min.] to select the next parameter .
⇒ The next digits of the MAC address are displayed.

7. At levels L72 to L75, repeat step 5 until the MAC address has been completely read out.

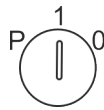

8. Set key switch to “1” to complete your input.
⇒ The MAC address has been read out.

7.6 Stocking the unit



CAUTION!

Danger to chilled goods due to contamination!

A dirty unit can contaminate the stored chilled goods. A contamination can damage or destroy the chilled goods.

- Do not store contaminated chilled goods in the unit.
- After any contamination, clean, disinfect and sterilise the unit and the chilled goods.
- Clean, disinfect and sterilise the unit in line with this instruction manual during operation.
- Comply with the applicable requirements for personal hygiene, disinfection and sterilisation.

During operation, the unit can be stocked with chilled goods at any time.

If the unit door is opened for more than 60 seconds during stocking, the door-open alarm is triggered & Chapter 10 “Alarms” on page 56.

Observe the following when stocking the unit:

- Make sure that the chilled goods match the intended purpose of the unit.
- Make sure that no infected or contaminated chilled goods are stored.
- Make sure that the chilled goods are stored in closed containers.
- Observe the maximum carrying capacity of the wire shelves and drawers.

Stocking the unit

- Before and during stocking, comply with the applicable regulations on personal hygiene.
- Before and during stocking, comply with the safety requirements relevant for the type of chilled goods.



CAUTION!

Danger of injury from broken drawers or shelves and falling chilled goods after overloading!

The carrying capacity of the shelves and drawers is limited. If the shelves and drawers are overloaded, they can break. There is a danger of cutting injuries caused by broken edges. Falling chilled goods can cause bruising.

- Load the drawers and shelves with a maximum of 100 kg/m² (as per DIN 58345/58371/58375).
- Store chilled goods in break-proof containers.
- Make sure that the chilled goods are stable and cannot fall over.

Requirements:

- The unit has been commissioned and the target temperature has been reached.
 - The cooling chain for the chilled goods as specified by the relevant manufacturer has been maintained.
1. ➤ Stock the unit with chilled goods.
 2. ➤ If the door is open for more than 60 seconds, deactivate the door-open alarm if necessary.

8 Cleaning and disinfection

Suitable disinfectants

The disinfectants listed in the following table have been tested by KIRSCH at the factory.

Adhere to the instruction manuals of the relevant manufacturers.

Tab. 8: Disinfectants

Disinfectants	Manufacturer
Incidin liquid	Ecolab Deutschland GmbH
Mikrozid AF liquid	Schülke & Mayr GmbH
Bacillol 30 Foam	Bode Chemie GmbH



Using other disinfectants

If disinfectants other than those mentioned above are used, test them at an inconspicuous location before their first use.

Use only acid-free disinfectants.

If in doubt, contact KIRSCH.

Cleaning and disinfecting the unit

Protective equipment: ■ Chemical-resistant protective gloves



CAUTION!

Danger to chilled goods due to contamination!

A dirty unit can contaminate the stored chilled goods. A contamination can damage or destroy the chilled goods.

- Do not store contaminated chilled goods in the unit.
- After any contamination, clean, disinfect and sterilise the unit and the chilled goods.
- Clean, disinfect and sterilise the unit in line with this instruction manual during operation.
- Comply with the applicable requirements for personal hygiene, disinfection and sterilisation.



CAUTION!

Danger of damage due to interrupted cooling chain during stock transfer!

If the cooling chain of the chilled goods is interrupted by a stock transfer, the prescribed storage conditions are no longer met. This can damage the chilled goods.

- Do not subject chilled goods to light during the stock transfer.
- Do not place chilled goods near radiators during the stock transfer.
- Make sure that chilled goods are stored in the replacement unit according to the specifications of the relevant manufacturer.

Requirements:

- The chilled goods have been transferred.
- The unit has been switched off ↪ *Chapter 7.2 “Switching off the unit” on page 41.*

1. ➤ Remove drawers and shelves.



NOTICE!

Danger of damage to the electrical system due to unsuitable cleaning agents!

The unit contains sensitive electronic parts. If the electronic parts come into contact with unsuitable cleaning agents, this can lead to a loss of function in the unit.

- Use cleaning agents free from sand and acid.
- Do not use chemical solvents.
- Do not bring cleaning agents into contact with electrical parts.

- 2.** ➤ Clean interior with warm water, dry it, and then disinfect it with a suitable disinfectant (↪ *“Disinfectants” on page 50*).
- 3.** ➤ Clean the drawers and shelves with washing-up liquid, dry them and disinfect them with a suitable disinfectant (↪ *“Disinfectants” on page 50*).
- 4.** ➤ Put drawers and shelves back in the unit.
- 5.** ➤ Wipe the door seal only with clear water and rub until thoroughly dry.
- 6.** ➤ Switch on the unit ↪ *Chapter 7.1 “Switching on the unit” on page 41.*

Cleaning the housing

Painted housing

- ▶ Treat the housing with cleaning and care products for painted surfaces.

Stainless steel housing

- ▶ Treat the housing with cleaning and care products for stainless steel.

9 Maintenance

Interval	Maintenance work
At least every 6 months	Check the condenser ↪ <i>Chapter 9.1 "Safety inspection" on page 53.</i>

9.1 Safety inspection



Perform safety inspection every two years!

The unit should be inspected by the owner in line with DGUV regulation 3 (formerly BGV A3).

Subject the unit to a safety inspection at regular intervals, but at least every two years.

Contents of the safety inspection

The safety inspection performed by the owner contains the following individual inspections and their documentation:

- 1 - Visual inspection
- 2 - Function check
- 3 - Temperature check
- 4 - Temperature alarm test
- 5 - Inspection of the condenser



NOTICE!

If there is any doubt that the unit is in the correct condition, immediately stop using the unit. To prevent unintended use, label the unit accordingly. Contact the service department ↪ *"Service contact" on page 4.*

Visual inspection

1. ➤ Check the entire unit for completeness, correct set-up and damage.
2. ➤ Check the following parts of the unit individually for damage:
 - Door handle
 - Interior
 - Door seal
3. ➤ If there is damage and functionality is not ensured, decommission the unit and contact the service department ↪ *"Service contact" on page 4.*

Function check


- Check that the following parts of the unit are functional:
 - Door handle

Checking the temperature


Special tool: ■ Calibrated temperature measuring device

1. ➤ Create a simulation of chilled goods for the sensor.
2. ➤ Attach measuring sensor in the interior at medium height.
3. ➤ Read temperature after 120 minutes.
4. ➤ After the check, clean and disinfect the unit ↪ *Chapter 8 "Cleaning and disinfection" on page 50.*

Testing the temperature alarm

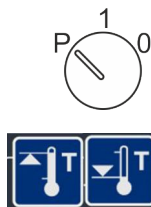

1. ➤ Turn unit key to position "P".
2. ➤ Press buttons [*Temperature warning max.*] and [*Temperature warning min.*] simultaneously and hold for approx. four seconds.
 - ⇒ The Display shows a flashing decimal point (.).
 - The test function starts, the electronic delay is switched off for 10 minutes.
3. ➤ Heat the monitoring sensor (↪ *Chapter 1.5.7 "Temperature sensor" on page 18*) (for example with your fingers).
4. ➤ Wait until the warning limit is exceeded and the buzzer sounds.
 - ⇒ The Display alternates between the current temperature and the error message.
5. ➤ Cool down the monitoring sensor (for example with a cooling spray).
6. ➤ Wait until the warning limit is exceeded and the buzzer sounds.
 - ⇒ The Display alternates between the current temperature and the two error messages (upper and lower temperature alarm).

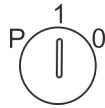

7. Turn unit key to position “1”.

⇒ The test function is completed, the electronic delay is switched on again.

The Display displays the current temperature of the chilled goods.



The test function is ended automatically after 10 minutes.

8. After the check, clean and disinfect the unit ↗ *Chapter 8 “Cleaning and disinfection” on page 50.*

Checking the condenser

Remove the dust from the condenser (Fig. 6) at least every six months in order to avoid affecting the performance of the cooling machine.

Requirement:

■ The back of the unit is accessible.

1. Dust off the condenser (for example with a brush or vacuum cleaner).

2. Check the condenser for visible damage and wear.

10 Alarms

10.1 Alarm functions

Alarm functions

If a function of the unit is faulty or defective, an alarm is triggered.

Every alarm is displayed visually as well as acoustically.

The display alternates between the visual alarm and the temperature. The message is displayed until the alarm is acknowledged. Acknowledging the alarm does not rectify the error.

The acoustic alarm is output as an alarm sound (hereafter referred to as “buzzer”).

The unit is equipped with the following alarm functions:






- Temperature alarm
- Door-open alarm
- Power failure warning
- Alarm in case of defective display and control unit

Depending on the cause that triggered the alarm, the following measures are required:

- 1 - Deactivate buzzer ↪ *“Deactivating the buzzer” on page 57*
- 2 - Acknowledge alarm ↪ *“Acknowledging the alarm” on page 57*
- 3 - Inform the service department ↪ *“Service contact” on page 4*

Tab. 9: Alarm functions (overview)

Alarm function	Display	Buzzer	Cause	Measure
Temperature alarm		✓	<ul style="list-style-type: none"> ■ The temperature is over the temperature warning limit. ■ The remote warning contact has been triggered. 	<ul style="list-style-type: none"> ■ Deactivate buzzer. ■ Determine cause and rectify.
		✓	<ul style="list-style-type: none"> ■ The temperature is below the temperature warning limit. ■ The remote warning contact has been triggered. 	
Door-open alarm		✓	<ul style="list-style-type: none"> ■ The door is open for more than 60 seconds. ■ The door-open alarm is not forwarded via the remote warning contact. 	<ul style="list-style-type: none"> ■ Deactivate buzzer. ■ Close door.

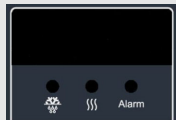Alarm function	Display	Buzzer	Cause	Measure
Power failure warning		✓	<ul style="list-style-type: none"> ■ The power supply of the unit has failed. ■ The remote warning contact has been triggered. 	<ul style="list-style-type: none"> ■ Deactivate buzzer. ■ Determine cause of power failure and rectify. ■ The monitoring unit remains in operation for approx. 72 hours. ■ Acknowledge alarm.
Alarm in case of defective battery		✓	<ul style="list-style-type: none"> ■ The power supply for the temperature documentation has failed. ■ The alarm function has failed. 	<ul style="list-style-type: none"> ■ Inform the service department. ■ Replace the battery  "Service contact" on page 4.
Alarm in case of defective display and control unit		✓	<ul style="list-style-type: none"> ■ The display and control unit is defective. 	<ul style="list-style-type: none"> ■ Unplug power plug. ■ Inform the service department  "Service contact" on page 4.

10.2 Handling alarms

Deactivating the buzzer

The buzzer sounds.



- Press button [Reset].
⇒ The buzzer is deactivated.
- Determine cause of the alarm and rectify.

Otherwise the buzzer will sound again every 30 minutes.

Acknowledging the alarm

The buzzer is deactivated, the alarm cause has been rectified, but the alarm continues to appear on the display.



- Press button [Reset].
⇒ The alarm has been acknowledged. The display shows the temperature of the chilled goods.

11 Status displays and error messages

11.1 Status displays

Status displays provide information to the user (for example regarding an ongoing defrosting process).

A status display is not accompanied by an acoustic signal (hereafter referred to as “buzzer”)

A status display does not require any immediate action by the user.

Tab. 10: Status displays

Display	Buzzer	Description	Measure	Unit key
	–	LED display [Defrosting] lights up: <ul style="list-style-type: none"> ■ The defrosting process is active. 	–	–
	–	LED display [Alarm] lights up: <ul style="list-style-type: none"> ■ One or more alarms have been triggered (collective alarm). 	–	–
	–	Standby display: <ul style="list-style-type: none"> ■ The unit is connected to the power grid and the key switch is set to position “0”. 	<ul style="list-style-type: none"> ■ Turn key switch to position “1” to switch on the unit. 	✓

11.2 Error messages

Do not make repairs yourself



WARNING!

Danger due to incorrect repairs or changes!

Incorrect repairs or changes can cause serious injury (e.g. electric shocks) and damage (e.g. fire, damage to chilled goods).

- Have repairs performed by the service department.
- Use KIRSCH replacement parts.
- Do not make independent additions or changes to the unit.
- If in doubt, contact KIRSCH.

Transferring chilled goods



NOTICE!

Danger to chilled goods due to defective or faulty unit!

A defect or fault in the unit means that its cooling performance is no longer ensured. Reduced cooling performance can cause considerable damage to chilled goods.

- Select an alternative storage location for the chilled goods.
- Ensure operating and storage conditions.
- Transfer chilled goods to new location.

Occurrence of error messages

Error messages indicate a malfunction of the unit.

Error messages and the temperature display alternate on the display.

If there is more than one error, the errors are shown one after the other on the display.

In addition to the information on the display, the buzzer sounds to report the error.

The unit indicates the following errors visually and acoustically:

- Unit errors
- Software errors

When error messages occur, proceed as described below:

Procedure in case of error messages

1. ➤ Deactivate buzzer.
2. ➤ Assess error display according to the table (↪ Chapter 11 “Status displays and error messages” on page 58).
3. ➤ Perform the recommended measures.
4. ➤ Acknowledge the alarm messages.



For repairs, contact the service department:

The following company has been commissioned and authorised by KIRSCH to provide service for the unit: ↪ “Service contact” on page 4



CAUTION!

Danger of damage due to interrupted cooling chain during stock transfer!

If the cooling chain of the chilled goods is interrupted by a stock transfer, the prescribed storage conditions are no longer met. This can damage the chilled goods.

- Do not subject chilled goods to light during the stock transfer.
- Do not place chilled goods near radiators during the stock transfer.
- Make sure that chilled goods are stored in the replacement unit according to the specifications of the relevant manufacturer.











Meaning of “X” for error and status messages

X is not shown on the display.


- Instead, the display shows a number that describes the relevant part.

Tab. 11: Error messages of the unit






Display	Buzzer	Description	Measure	Unit key required
	✓	Sensor X: <ul style="list-style-type: none"> ■ Error or short circuit in the relevant sensor. ■ The cryostat is running in the emergency program. 	<ul style="list-style-type: none"> ■ Contact the service department. 	-
	✓	Sensor X: <ul style="list-style-type: none"> ■ Error or break in the relevant sensor. ■ Cryostat is running in the emergency program. 	<ul style="list-style-type: none"> ■ Contact the service department. 	-
	✓	Fan X: <ul style="list-style-type: none"> ■ Speed of the relevant fan too low. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
	✓	Fan X: <ul style="list-style-type: none"> ■ Speed of the relevant fan too high. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
	✓	Fan: <ul style="list-style-type: none"> ■ Fan does not reach required minimum speed after an unit restart. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
	✓	Fan: <ul style="list-style-type: none"> ■ Difference between the speeds of the fans is too great. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
	✓	Relay X: <ul style="list-style-type: none"> ■ Defect in the relevant relay. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
	✓	Relay X: <ul style="list-style-type: none"> ■ Defect in the relevant relay. ■ The temperature of the chilled goods can fluctuate. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-

Error messages

Display	Buzzer	Description	Measure	Unit key required
PAR	✓	<p>Synchronisation error:</p> <ul style="list-style-type: none"> ■ Synchronisation error between control unit and monitoring circuit. ■ No secured function of the cooling controller. 	<ul style="list-style-type: none"> ■ Set key switch to "0". ■ Unplug power plug and switch on again ↪ <i>Chapter 7.1 "Switching on the unit" on page 41.</i> 	✓
Con	✓	<p>Connection problem:</p> <ul style="list-style-type: none"> ■ Synchronisation error between control unit and monitoring circuit. ■ No secured function of the cooling controller. 	<ul style="list-style-type: none"> ■ Transfer chilled goods to new location. ■ Contact the service department. 	-
IEC	✓	<p>Control error:</p> <ul style="list-style-type: none"> ■ Error during self-test in the cooling controller. ■ The monitoring circuit takes over temperature control. 	<ul style="list-style-type: none"> ■ Contact the service department. 	-
door	✓	<p>Door-open alarm:</p> <ul style="list-style-type: none"> ■ Door is open for more than 60 seconds. 	<ul style="list-style-type: none"> ■ Close door. 	-
bATT	✓	<p>Battery defective:</p> <ul style="list-style-type: none"> ■ Battery must be replaced. ■ The temperature documentation and alarm in the event of power failure fail. 	<ul style="list-style-type: none"> ■ Switch off alarm. ■ Inform the service department. ■ Replace the battery ↪ <i>"Service contact" on page 4.</i> 	-
PF	✓	<p>Power failure:</p> <ul style="list-style-type: none"> ■ Unit is not cooling. ■ Alarm is active. ■ The remote warning contact is triggered. 	<ul style="list-style-type: none"> ■ Check power supply. 	-
EXI	✓	<p>Temperature alarm (high):</p> <ul style="list-style-type: none"> ■ The upper temperature warning limit is reached or exceeded (for example due to very warm chilled goods or because the door was open too long). 	<ul style="list-style-type: none"> ■ View and check the temperature warning limit ↪ <i>Chapter 6.2.2.2 "Displaying and changing the temperature warning limits" on page 38.</i> ■ If necessary, correct values of the temperature warning limits. ■ Monitor the temperature progression. ■ If it does not normalise, contact the service department. ■ Transfer chilled goods to new location. 	✓

Display	Buzzer	Description	Measure	Unit key required
	✓	<p>Temperature alarm (low)</p> <ul style="list-style-type: none"> ■ The lower temperature warning limit is reached or exceeded (for example after the door has been opened for a longer period while the cooling machine cools down the unit). ■ Safety equipment triggers, monitoring circuit has switched off the cooling machine. 	<ul style="list-style-type: none"> ■ View and check the temperature warning limit ↪ <i>Chapter 6.2.2.2 “Displaying and changing the temperature warning limits” on page 38.</i> ■ If necessary, correct values of the temperature warning limits. ■ Monitor the temperature progression. ■ If it does not normalise, contact the service department. ■ Transfer chilled goods to new location. 	✓

Tab. 12: Error messages and status displays of the PC-KIT-STICK

Display	Buzzer	Description	Measure	Unit key
	–	<p>Status display:</p> <ul style="list-style-type: none"> ■ Copying process is running. ■ XX represents the progress of the copying process in percent. 	<ul style="list-style-type: none"> ■ Do not remove USB flash drive. 	–
	–	<p>Status display:</p> <ul style="list-style-type: none"> ■ Copying process is complete. 	<ul style="list-style-type: none"> ■ Remove USB flash drive. 	–
	–	<p>Status display:</p> <ul style="list-style-type: none"> ■ Memory of the USB flash drive is full. 	<ul style="list-style-type: none"> ■ Use USB flash drive with sufficient memory capacity. 	–
	–	<p>Error message:</p> <ul style="list-style-type: none"> ■ The USB flash drive cannot be read. 	<ul style="list-style-type: none"> ■ Use a USB flash drive formatted in FAT-32 format. 	–
	–	<p>Error message:</p> <ul style="list-style-type: none"> ■ USB flash drive was removed while copying process was ongoing. 	<ul style="list-style-type: none"> ■ Remove USB flash drive and try again after 1 minute. 	–

12 Decommissioning and disposal

12.1 Decommissioning unit

Decommissioning

1. ▶ Switch off unit.
2. ▶ Transfer chilled goods to new location.
3. ▶ Unplug power plug.
4. ▶ Cut through connecting cable.
5. ▶ Remove or destroy locks.
6. ▶ Remove door.

12.2 Disposing of the unit

Disposing of the battery



ENVIRONMENT!

Danger to the environment due to incorrect disposal of the battery!

If the battery is disposed of separately, proceed as follows:

- Do not damage, burn or short-circuit the battery.
- Dispose of the battery in accordance with regional regulations.
- If in doubt, consult the local communal authority or special disposal specialists regarding environmentally friendly disposal.

Disposing of the unit



ENVIRONMENT!

Danger to the environment due to incorrect disposal of the unit!

If substances hazardous to the environment are handled incorrectly, and especially if they are disposed of incorrectly, this can cause serious damage to the environment.

- Do not independently remove and dispose of the cooling machine.
- If substances hazardous to the environment (refrigerant, for example) accidentally enter the environment, take suitable measures immediately. If in doubt, contact the responsible communal authority to report the damage and enquire about the measures to be taken.
- Dispose of the unit in accordance with the regional regulations for electrical and electronic units.
- If in doubt, consult the local communal authority or special disposal specialists regarding environmentally friendly disposal.

13 Appendix

Depending on the model, the appendix contains the following applicable documents:

- Declaration of conformity
- Technical data
- Installation drawing

13.1 Declaration of conformity

EC Declaration of Conformity

We,

Philipp Kirsch GmbH
Im Lossenfeld 14
77731 Willstätt-Sand
Germany

declare that the devices described below comply with the protection requirements of the directives and standards below at the time that the devices were placed on the market.

Manufacturer	Device category	Type	Serial numbers from serial number onwards
KIRSCH	Refrigerator	LABO 100 PRO-ACTIVE	100 33 25000
KIRSCH	Refrigerator	LABO 100 PRO-ACTIVE	100 05 25000
KIRSCH	Refrigerator	LABO 126 PRO-ACTIVE	125 05 25000
KIRSCH	Refrigerator	LABO 126 PRO-ACTIVE	125 34 25000
KIRSCH	Refrigerator	LABO 288 PRO-ACTIVE	280 34 25000
KIRSCH	Refrigerator	LABO 288 PRO-ACTIVE	280 15 25000
KIRSCH	Refrigerator	LABO 340 PRO-ACTIVE	330 31 25000
KIRSCH	Refrigerator	LABO 340 PRO-ACTIVE	330 12 25000
KIRSCH	Refrigerator	LABO 468 PRO-ACTIVE	460 31 25000
KIRSCH	Refrigerator	LABO 468 PRO-ACTIVE	460 05 25000
KIRSCH	Refrigerator	LABO 520 PRO-ACTIVE	500 30 25000
KIRSCH	Refrigerator	LABO 520 PRO-ACTIVE	500 01 25000
KIRSCH	Refrigerator	LABO 520 CR PRO-ACTIVE	500 42 25000
KIRSCH	Refrigerator	LABO 520 CR PRO-ACTIVE	500 02 25000
KIRSCH	Refrigerator	LABO 720 PRO-ACTIVE	700 30 25000
KIRSCH	Refrigerator	LABO 720 PRO-ACTIVE	700 01 25000
KIRSCH	Refrigerator	LABO 720 CR PRO-ACTIVE	700 36 25000
KIRSCH	Refrigerator	LABO 720 CR PRO-ACTIVE	700 02 25000
KIRSCH	Refrigerator	LABO 720 CHROMAT PRO-ACTIVE	700 32 25000
KIRSCH	Refrigerator	LABO 720 CHROMAT PRO-ACTIVE	700 20 25000
KIRSCH	Refrigerator	LABO 288 ULTIMATE	280 70 25000
KIRSCH	Refrigerator	LABO 288 ULTIMATE	280 80 25000
KIRSCH	Refrigerator	LABO 340 ULTIMATE	340 70 25000
KIRSCH	Refrigerator	LABO 340 ULTIMATE	340 80 25000
KIRSCH	Refrigerator	LABO 468 ULTIMATE	460 70 25000
KIRSCH	Refrigerator	LABO 468 ULTIMATE	460 80 25000
KIRSCH	Refrigerator	LABO 520 ULTIMATE	500 70 25000
KIRSCH	Refrigerator	LABO 520 ULTIMATE	500 80 25000
KIRSCH	Refrigerator	LABO 720 ULTIMATE	700 70 25000
KIRSCH	Refrigerator	LABO 720 ULTIMATE	700 80 25000
KIRSCH	Refrigerator	LABO 720 CHROMAT ULTIMATE	700 74 25000
KIRSCH	Refrigerator	LABO 720 CHROMAT ULTIMATE	700 87 25000

Directives

- RoHS-Directives 2011/65/EG
- EC Low Voltage Directives 2014/35/EU
- Electromagnetic Compatibility 2014/30/EU
- Machinery Directives 2006/42/EG

Standards

- DIN EN ISO 9001:2015
- DIN 13221:2016-08

Harmonised standards

- EN 61010-1:2010+A1:2019
- EN 61010-2-11:2017
- EN 60601-1-2:2015
- EN/IEC 61000-3-2:2014
- EN/IEC 61000-3-3:2013

Declaration of conformity

We also declare that the devices described below comply with the protection requirements of the directives and standards below at the time that the devices were placed on the market.

Manufacturer	Device category	Type	Serial numbers from serial number onwards
KIRSCH	Freezer	FROSTER LABO 330 PRO-ACTIVE	320 01 25000
KIRSCH	Freezer	FROSTER LABO 530 PRO-ACTIVE	500 08 25000
KIRSCH	Freezer	FROSTER LABO 530 CR PRO-ACTIVE	500 09 25000
KIRSCH	Freezer	FROSTER LABO 730 PRO-ACTIVE	700 04 25000
KIRSCH	Freezer	FROSTER LABO 730 CR PRO-ACTIVE	700 05 25000
KIRSCH	Freezer	FROSTER LABO 330 ULTIMATE	320 80 25000
KIRSCH	Freezer	FROSTER LABO 530 ULTIMATE	500 84 25000
KIRSCH	Freezer	FROSTER LABO 730 ULTIMATE	700 84 25000

Directives

- RoHS-Directives 2011/65/EG
- EC Low Voltage Directives 2014/35/EU
- Electromagnetic Compatibility 2014/30/EU
- Machinery Directives 2006/42/EG

Standards

- DIN EN ISO 9001:2015

Harmonised standards

- EN 61010-1:2010+A1:2019
- EN 61010-2-11:2017
- EN 60601-1-2:2015
- EN/IEC 61000-3-2:2014
- EN/IEC 61000-3-3:2013



Willstätt, 05.05.2021
 Dr. Jochen Kopitzke
 Managing Director




13.2 Technical data

	LABO 100 PRO-ACTIVE	LABO 126 PRO-ACTIVE	LABO 288 PRO-ACTIVE/ ULTIMATE*	LABO 340 PRO-ACTIVE/ ULTIMATE*
Capacity in litres	95	95	280	330
Temperature setting approx. in °C	approx. 0 to +15	approx. 0 to +15	approx. 0 to +15	approx. 0 to +15
Voltage	220–240	220–240	220–240	220–240
Frequency in Hz	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50
Refrigerant quantity in grams	Frigen = 100 Eco = 42	Frigen = 70 Eco = 32	Frigen = 80 Eco = 30	Frigen = 160 Eco = 40
Power consumption in watts	Frigen = 160 Eco = 76	Frigen = 160 Eco = 105	Frigen = 155 Eco = 88	Frigen = 150 Eco = 110
Normal consumption in kWh/24	Frigen = 0.86 Eco = 0.64	Frigen = 0.86 Eco = 0.58	Frigen = 1.00 Eco = 0.68	Frigen = 1.20 Eco = 0.65
Perm. ambient tem- perature in °C	+10 to +38	+10 to +38	+10 to +38	+10 to +38
Exterior dimensions including wall dis- tance (WxDxH) in cm	54 x 54 x 82	54 x 53.5 x 88.8	67 x 75 x 124/ 67 x 75 x 131*	67 x 69 x 181 67 x 69 x 188*
Interior dimensions (WxDxH) in cm	44 x 43 x 47 (usable depth 5 cm, usable width 2 cm less)	44 x 42 x 67 (usable depth at top 5 cm, at bottom 15 cm less, usable width 2 cm less)	53 x 50 x 100 (usable depth at top 5 cm, at bottom 13 cm less)	53 x 45 x 128 (usable depth 5 cm less)
Exterior dimensions with door open 90° (WxD) in cm	54 x 106	54 x 106	67 x 130	67 x 127
Shelf size (WxD) in cm	43 x 33	43 x 33	52 x 39	53 x 39
Clear drawer dimen- sion (WxDxH) in cm	40 x 32 x 5.6	-	-	50 x 32 x 5.6
Max. load capacity drawer/shelf in kg	13 / 25	13 / 25	16 / 40	16 / 40
Net/gross weight in kg	46/50	47 / 57	78 / 85	88 / 100
Noise emission (for Eco! R600a) in dB(A)	41	40	41	39

	LABO 468 PRO-ACTIVE/ ULTIMATE*	LABO 520 PRO-ACTIVE/ ULTIMATE*	LABO 720 PRO-ACTIVE/ ULTIMATE*	LABO 720 CHROMAT PRO-ACTIVE/ ULTIMATE*
Capacity in litres	460	500	700	700
Temperature setting approx. in °C	approx. 0 to +15	approx. 0 to +15	approx. 0 to +15	approx. 4 to +15
Voltage	220–240	220–240	220–240	220–240
Frequency in Hz	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50	Frigen = 50/60 Eco = 50
Refrigerant quantity in grams	Frigen = 180 Eco = 65	Frigen = 350 Eco = 90	Frigen = 350 Eco = 90	Frigen = 350 Eco = 90
Power consumption in watts	Frigen = 370 Eco = 234	Frigen = 370 Eco = 250	Frigen = 370 Eco = 250	Frigen = 370 Eco = 250
Normal consumption in kWh/24	Frigen = 1.40 Eco = 1.00	Frigen = 1.40 Eco = 1.40	Frigen = 1.46 Eco = 1.46	Frigen = X.XX Eco = 1.63
Perm. ambient tem- perature in °C	+10 to +38	+10 to +38	+10 to +38	+10 to +38
Exterior dimensions including wall dis- tance (WxDxH) in cm	74 x 80 x 181/ 74 x 80 x 188*	77 x 76 x 195.5	77 x 98 x 195.5	77 x 98 x 195.5
Interior dimensions (WxDxH) in cm	60 x 57 x 128 (usable depth 5 cm less)	62 x 57 x 140 (usable depth 11 cm, usable height 13 cm, usable width 2 cm less)	62 x 77 x 140 (usable depth 11 cm, usable height 13 cm, usable width 2 cm less)	62 x 77 x 140 (usable depth 11 cm, usable height 13 cm, usable width 2 cm less)
Exterior dimensions with door open 90° (WxD) in cm	74 x 144	74 x 144	77 x 165	77 x 165
Shelf size (WxD) in cm	59 x 50	59 x 45	59 x 65	59 x 65
Clear drawer dimen- sion (WxDxH) in cm	-	-	-	-
Max. load capacity drawer/shelf in kg	24 / 40	- / 40	- / 40	- / 40
Net/gross weight in kg	109 / 122	120 / 148	145 / 174	145 / 174
Noise emission (for Ecool R600a) in dB(A)	41	42	42	42

	FR LABO 330 PRO-ACTIVE/ ULTIMATE*	FR LABO 530 PRO-ACTIVE/ ULTIMATE*	FR LABO 730 PRO-ACTIVE/ ULTIMATE*
Capacity in litres	300	500	700
Temperature setting approx. in °C	approx. -10 to -30	approx. -10 to -30	approx. -10 to -30
Voltage	220–240	220–240	220–240
Frequency in Hz	Frigen = 50/60	Frigen = 50/60	Frigen = 50/60
Refrigerant quantity in grams	Frigen = 440	Frigen = 380	Frigen = 380
Power consumption in watts	Frigen = 590	Frigen = 605	Frigen = 605
Normal consumption in kWh/24	Frigen = 3.69	Frigen = 7.22	Frigen = 7.51
Perm. ambient tem- perature in °C	+10 to +32	+10 to +32	+10 to +32
Exterior dimensions including wall dis- tance (WxDxH) in cm	74 x 77 x 159/ 74 x 77 x 166*	77 x 76 x 195.5	77 x 97 x 195.5
Interior dimensions (WxDxH) in cm	59.5 x 53 x 95 (usable depth 10 cm, usable height 13 cm, usable width 2 cm less)	60 x 57 x 140 (usable depth 11 cm, usable height 13 cm, usable width 2 cm less)	62 x 78 x 140 (usable depth 11 cm, usable height 13 cm, usable width 2 cm less)
Exterior dimensions with door open 90° (WxD) in cm	74 x 141	77 x 144	77 x 165
Shelf size (WxD) in cm	57 x 42	59 x 65	59 x 65
Clear drawer dimen- sion (WxDxH) in cm	-	-	-
Max. load capacity drawer/shelf in kg	40	40	40
Net/gross weight in kg	127 / 142	140 / 183	165 / 200
Noise emission (for Ecool R600a) in dB(A)	51.2	53	53

13.3 Installation drawing

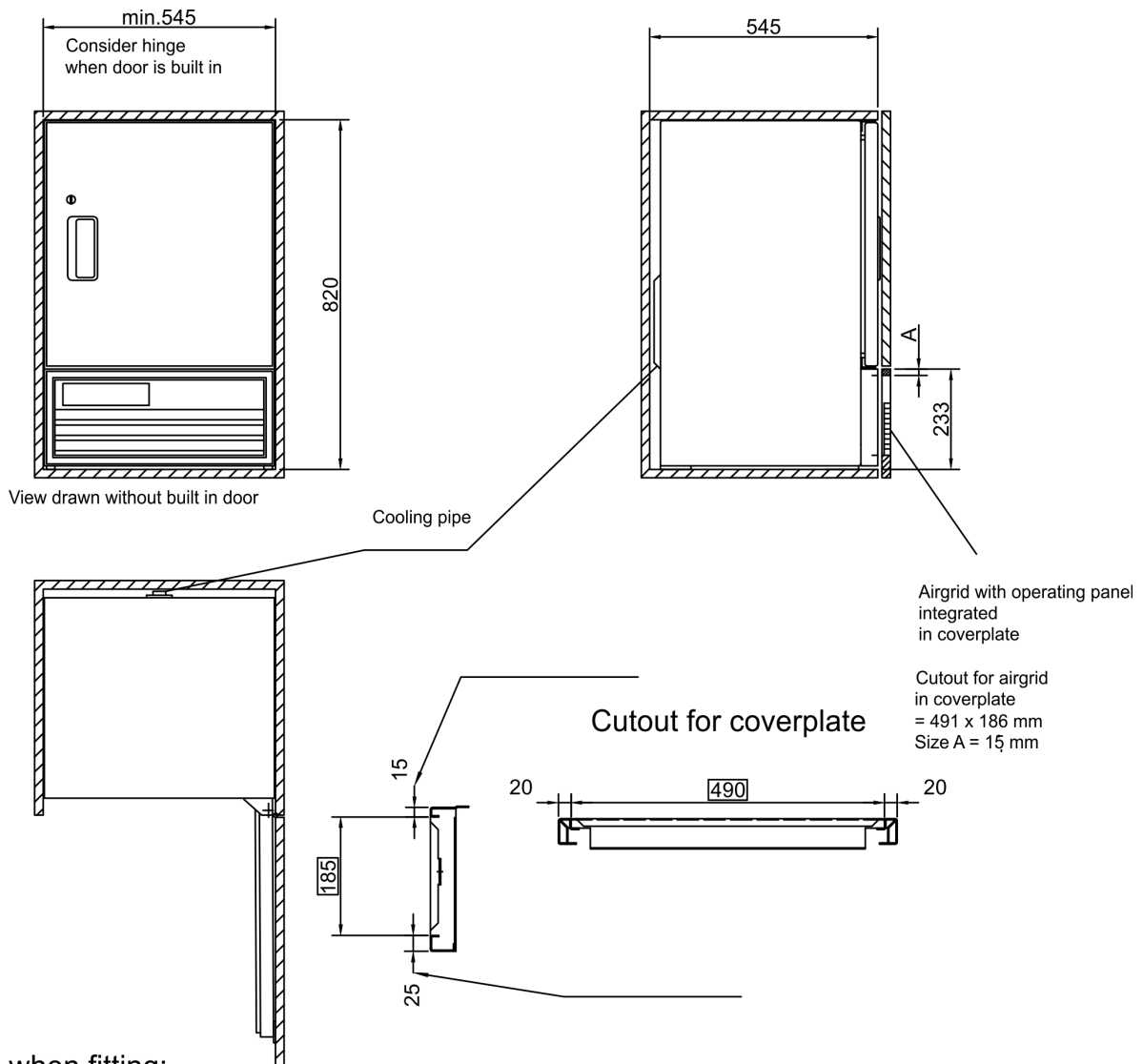


Ensure when fitting:


Supply and exhaust air air takes place via front-sided airgrid.
For that do not block with subjects or even blind it,
so that the ventilation of the cooling machine keeps warranted.

Drawing.Nr.: 225-033-1

Fig. 10: Installation drawing for LABO 100 PRO-ACTIVE

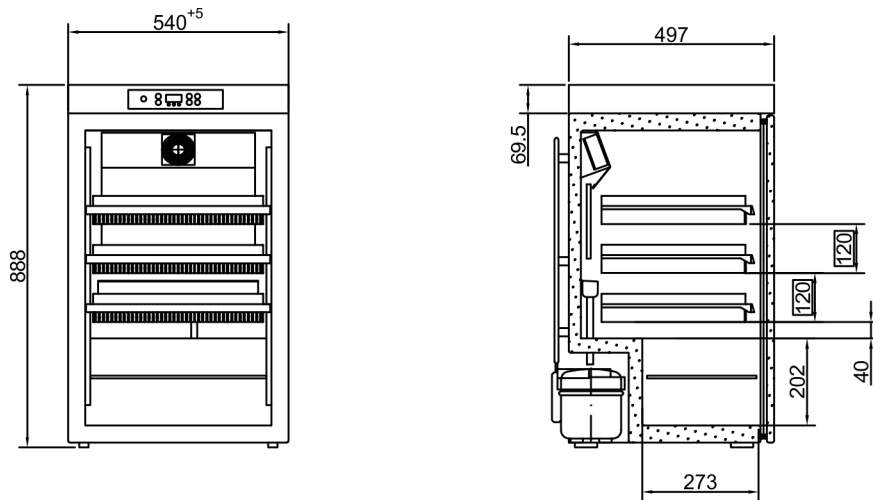

Ansicht ohne Türe gezeichnet

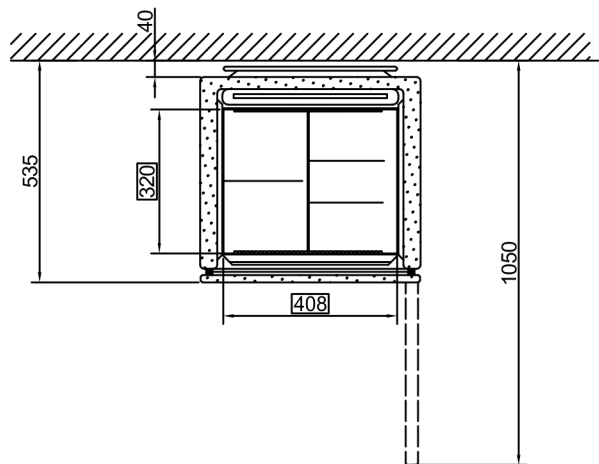

Bei Einbau beachten:

Luftdurchgang oben und unten :

ohne Kühlmaschinenventilator: mind. 250 cm²
 mit Kühlmaschinenventilator: mind. 80 cm²

Die mit umrandeten Maße
 geben die nutzbare Breite / Tiefe / Höhe an.

Fig. 11: Installation drawing for LABO 126 PRO-ACTIVE

14 Index

A	
About this instruction manual	3
Additional defrosting	16
Alarm functions	
Alarms (overview)	56
Appendix	66
B	
Button functions (overview)	12
C	
Change of location	29
Character	23
Checking the condenser	55
Circulation cooling	15
Cleaning	
Housing	52
Interior	50
Commissioning	
Process	34
Programming	34
Staff qualification	34
Condenser	16
Connection	33
Contact	4
Cooling	15
Cooling machine	15
Customer service	4
D	
Declaration of conformity	67
Decommissioning	64
Defrosting (automatic)	16
Defrosting (manual)	16
Disinfectants	50
Disinfection	50
Display	12
Disposal	64
Door key	14
Drawers	14
E	
Error messages (overview)	60
F	
Final decommissioning	31
Foreseeable misuse	23
Function check	54
H	
Humidity	40
I	
Images on the display (overview)	12
Installation	32
Installation drawing	72
Interfaces	14
Interior fittings	14
K	
Key switch	12
L	
LAN interface	14
M	
Manufacturer's address	4
O	
Operation	12
Retrieving/erasing the temperature memory	42
Setting up PC-KIT-NET	43
SuperFrost function	43
Switching off	41
Switching on	41
Other applicable documents	66
P	
PC-KIT-NET	14, 43
Potential-free contact	14
Protective equipment	28

Protective gloves	28	TCP/IP module	14, 43
Purpose	23	Technical data	69
R		Temperature alarm test	54
Recommissioning	31	Temperature check	54
Residual risks	24	Temperature display	17
S		Temperature memory	17
Safety	22	Temperature sensor (overview)	18
Safety inspection	53	Temperature warning limits	
Service department	4	(overview)	38
Set-up	32	changing	39
Set-up conditions	32	displaying	38
Shelves	14	Transport	29
Staff	27	Type plate	26
Status displays (overview)	58	U	
Storage	3	Unit door	14
Switching off	41	Glass door (optional)	20
Switching on	41	Unit key	14
Symbols	22	USB port	12, 14
T		V	
Target temperature	35	Visual inspection	53