



HandyStep® S

The uncomplicated repetitive pipette
for complicated liquids

BRAND. For lab. For life.®

- + High viscosity, high density, high vapor pressure – no problem, thanks to the positive displacement principle
- + Compatible with BRAND PD-Tips and other positive displacement tips
- + Robust, stable and functional



Get to know the HandyStep® S

The HandyStep® S repetitive pipette combines the application advantages of the positive displacement principle with easy operation, robust design and an attractive price-to-performance ratio.

Thanks to the positive displacement principle, the HandyStep® S is often the easy solution for liquids that are difficult to pipette. This includes highly viscous or wetting media, high density media, foaming media such as detergents or highly volatile liquids.

Its balanced weight distribution and the practical ergonomic positioning of operating elements makes work more efficient, even when pipetting long series. The ergonomic finger rest

relieves strain on the hand, the positioning of the dispensing lever is easy on the thumb, and the stroke can be set with one hand, even while wearing gloves. Tips can be ejected with the press of a button, without having to touch the tips, making tip ejection contamination-free.

The HandyStep® S repetitive pipette works with the original BRAND PD-Tips II. Tips from many other manufacturers can also be used, providing maximum flexibility.

The simple, robust design and purely mechanical, battery-free operation make the HandyStep® S a reliable lab partner. Discover the application possibilities of the HandyStep® S.

HandyStep® S



- + Volume range from 2 µl to 5 ml
- + Adjustable to 59 volumes and up to 49 dispensing steps
- + Weighs just 108 g
- + Easy tip connection – simply insert the PD-Tip

At a glance: Advantages of the HandyStep® S

The HandyStep® S is a robust repetitive pipette that uses the positive displacement principle. With the PD-Tips, the HandyStep® S covers a wide volume range from 2 µl to 5 ml

and can be used in many different applications, including diagnostics, molecular biology and environmental analysis.

Precise

Clear, large stroke selection indicator – for fast, simple checking.

Perfectly placed

Optimally positioned dispensing lever, easy on the thumb.

Safe

At the push of a button! Ergonomic eject button for contact-free, safe ejection of contaminated dispenser tips.

Straightforward

Double-sided volume table on the back side for quick adjustments

Rapid

One-handed volume adjustment for right- and left-handed operators – even with gloves

Comfortable


Ergonomic finger loop to keep the hand relaxed, important when conducting long series

Simplified

Just one lever for locking the piston and filling the PD-Tip //

Versatile

In addition to PD-Tips from BRAND, compatible dispenser tips from other manufacturers can also be used.

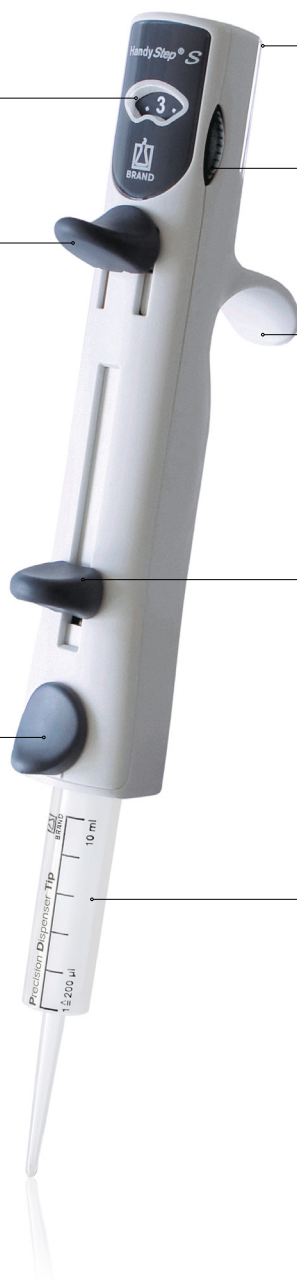

Volume table HandyStep® S with PD-Tips //

Offers 59 different partial volumes with different numbers of dispensing steps, depending on the PD-Tip size and stroke setting number used.

Setting	Tip size (ml)										Steps
	0.1	0.5	1	1.25	2.5	5	10	12.5	25	50	
1	2	10	20	25	50	100	200	250	500	1000	49
1.5	3	15	30	37.5	75	150	300	375	750	1500	32
2	4	20	40	50	100	200	400	500	1000	2000	24
2.5	5	25	50	62.5	125	250	500	625	1250	2500	19
3	6	30	60	75	150	300	600	750	1500	3000	15
3.5	7	35	70	87.5	175	350	700	875	1750	3500	13
4	8	40	80	100	200	400	800	1000	2000	4000	11
4.5	9	45	90	112.5	225	450	900	1125	2250	4500	10
5	10	50	100	125	250	500	1000	1250	2500	5000	9

Volume (µl)



Repetitive dispensing



Dispensing small volumes



The rack mount also fits the Transferpette® S bench-top rack.



Fast volume selection with the volume table on the back side

Ordering information

HandyStep® S

repetitive pipette, DE-M marked


Items supplied:

performance certificate, operating manual, shelf/rack mount.



Pack of	Cat. No.
1	705110

Accessories for simple and efficient work

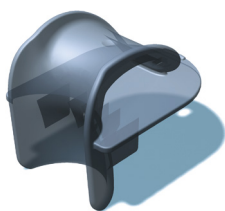

Shelf/rack mount

for 1 HandyStep® S repetitive pipette

Pack of	Cat. No.
1	705130



Wall mount

for 1 HandyStep® S repetitive pipette

Pack of	Cat. No.
1	705129

PD-Tips //

PD-Tips //, Precision dispenser tips with type code

Capacity ml	Pack of	non-sterile Cat. No.	BIO-CERT® LIQUID HANDLING STERILE Cat. No.
0.1	100	705700	705730
0.5	100	705702	705732
1	100	705704	705734
1.25	100	705706	705736
2.5	100	705708	705738
5	100	705710	705740
10	100	705712	705742
12.5	100	705714	705744
25	50 + 1 adapter	705716	705746*
50	25 + 1 adapter	705718	705748
Set	20 pieces per size: 0.5 1 1.25 2.5 5 10 12.5 ml	705720	



Adapter

for 25 and 50 ml PD-Tips //, PP,
autoclavable

Description	Pack of	Cat. No.
non-sterile	10	702398
BIO-CERT® LIQUID HANDLING STERILE	5	702399

* BIO-CERT®: pack of 25 + 1 adapter



Accuracy Table

The measurements were carried out with BRAND PD-Tips //. Distilled water was used as the test liquid. The nominal volume is the maximum volume printed on the PD-Tip. The tolerances specified in ISO 8655 are not exceeded.

PD-Tip //	Volume range	Accuracy A $\leq \pm$ %			Coefficient of variation CV \leq %		
		Stroke setting \pm % of nominal volume1 \pm 1 \pm 2% 3 \pm 6% 5 \pm 10%	Stroke setting \pm % of nominal volume1 \pm 1 \pm 2% 3 \pm 6% 5 \pm 10%	Stroke setting \pm % of nominal volume1 \pm 1 \pm 2% 3 \pm 6% 5 \pm 10%			
0.1 ml	2 μ l - 10 μ l	4.0	2.4	1.6	6.0	3.0	2.0
0.5 ml	10 μ l - 50 μ l	2.5	1.5	1.0	2.5	1.5	1.0
1 ml	20 μ l - 100 μ l	2.5	1.5	1.0	2.0	1.2	0.8
1.25 ml	25 μ l - 125 μ l	2.5	1.4	0.9	2.0	1.1	0.7
2.5 ml	50 μ l - 250 μ l	1.8	1.1	0.7	1.5	0.9	0.6
5 ml	100 μ l - 500 μ l	1.8	1.1	0.7	1.5	0.9	0.7
10 ml	200 μ l - 1000 μ l	1.8	1.1	0.7	2.0	1.2	0.8
12.5 ml	250 μ l - 1250 μ l	1.8	1.1	0.8	3.2	2.0	1.4
25 ml	500 μ l - 2500 μ l	1.5	0.9	0.6	3.0	1.5	1.0
50 ml	1000 μ l - 5000 μ l	1.5	0.8	0.5	5.0	1.8	1.2

BRAND GMBH + CO KG

P.O. Box 1155 | 97861 Wertheim | Germany

T +49 9342 808 0 | F +49 9342 808 98000 | info@brand.de | www.brand.de



BRAND. For lab. For life.®

HandyStep® S is made in Germany under BRAND patents US 8,652,4 18B2; EP 2 428 271 B1 and EP 2 415 525 B1. PD-Tips // possess features according to the BRAND patents USD825,750S and CN304542771S as well as the European Design Patents EM003763564 and EM004045847.

BIO-CERT®, BRAND®, BRAND. For lab. For life.®, HandyStep®, as well as the BRAND figurative mark are registered trademarks or trademarks of BRAND GMBH + CO KG, Germany. All other trademarks mentioned or depicted here are the property of the respective owners.

Our technical literature is intended to inform and advise our customers. However, the validity of general empirical values, and of results obtained under test conditions, for specific applications depends on many factors beyond our control. Please appreciate, therefore, that no claims can be derived from our advice. The user is responsible for checking the appropriateness of the product for any particular application.

California Residents: For more information concerning California Proposition 65, please refer to www.brand.de/calprop65.

Subject to technical modification without notice. Errors excepted.

994331 © 2022 BRAND GMBH + CO KG | Printed in Germany | 0422



Find accessories and replacement parts, user manuals, test instructions (SOP) and product videos at shop.brand.de



Further information on products and applications can be found on our YouTube channel: [mylabBRAND](https://www.youtube.com/mylabBRAND)



 **BRANDGROUP**

BRAND (Shanghai) Trading Co., Ltd.
Shanghai, China

Tel.: +86 21 6422 2318
info@brand.com.cn
www.brand.cn.com

BRAND Scientific Equipment Pvt. Ltd.
Mumbai, India

Tel.: +91 22 42957790
customersupport@brand.co.in
www.brand.co.in

BrandTech® Scientific, Inc.
Essex, CT. United States of America

Tel.: +1 860 767 2562
info@brandtech.com
www.brandtech.com