



LATERAL FLOW TEST KIT

for the detection of milk proteins in food, cip solutions and working surfaces

ProGnosis Biotech S.A. is ISO 9001:2015 certified by TÜV Hellas (TÜV NORD).

Use only the current version of Product Data Sheet enclosed with the kit.

Total Milk Free test, E2810/E2830, is a lateral flow test that detects milk proteins in food products, cip solutions and working surfaces. The lateral flow kit contains all reagents required for the immunoassay method.

Matrices:

Bakery products, beverages, bread mix, cake mix, chocolate, cooked biscuit, dressings, ice cream, sauces, sausage, sorbets, soup mix

- Sample preparation: extraction/swab sampling
- Test time (incubation time after samples and reagents preparation): 5 min
- Calibration Range: 4-100ppm
: 1-32 μ g/100cm² (for surfaces)
- Shelf life: 12 months
- Storage: 4-30°C

***Note:** The same extract can, also, be used in all the products of **Allergen Free** series. If you want to check your sample for more than one allergen at the same time, you can use the same extract to test as many of the products you have available.*

Method characteristics

- The Limit of Detection (LOD) of the method is 2.0 ppm (2.0mg/kg) milk protein or 62.37 ppm milk or 6.237 ppm milk powder in food samples and CIP solutions.
- The Limit of Quantification (LOQ) of the method is 4.0 ppm (4.0mg/kg) milk protein or 124.75 ppm milk or 12.475 ppm milk powder in food samples and CIP solutions.
- **Surfaces** LOD: 0.5µg/100cm² milk protein on working surfaces. LOD was calculated based on our reference materials.
- **Surfaces** LOQ: 1.0µg/100cm² milk protein on working surfaces. LOQ was calculated based on our reference materials.
- The antibodies specifically detect caseins as well as β-lactoglobulin of cow's, sheep's, goat's and buffalo's milk.

1. Description

Total Milk Free test is a Lateral Flow test for the detection of milk proteins in food products, specially for those labeled as milk-free, CIP solutions and working surfaces.

2. General Information

Cow's milk contains 3.2 % proteins which consist of 10 % β-lactoglobulin (most abundant protein of whey) and 80 % caseins. Consumption of milk might be harmful for people who are allergic to milk. The most common protein that may cause an immune reaction to children is β- lactoglobulin whereas adults usually may have an allergic reaction due to caseins. The allergen can be present as an ingredient or as a contamination in raw and cooked products. Consumption of milk-containing food from allergic people might cause a broad range of symptoms, such as hives, itching, mild oral allergy or/and anaphylactic shock. According to the regulation (EU) No. 1169/2011 Annex II, milk is included in the list of allergens established by the European Food Safety Authority, and its presence must be indicated on the label. Similar regulations exist e.g. in the USA, Canada, Australia and New Zealand.

3. Principle of the method

The presence of milk in a sample is determined by the immunological detection of milk proteins. Antibodies specific to milk proteins are coated on the test line region (Test line) of the nitrocellulose membrane. During testing, antigens in the specimen react with the antibodies that are coated onto gold nanoparticles. The mixture migrates up the membrane to react with the antibodies immobilized on the membrane and generate a colored line in the test region T. The presence of the colored Test line indicates a positive result. In case of samples with a very high allergen concentration, the Test line fades or may not appear, giving us reduced or false negative result (hook effect phenomenon). For this purpose, a second line has been created (Hook line), whose intensity decreases as the amount of antigen increases and at very high concentrations it disappears either together with the Test line or before it. To serve as a procedural control, a colored line will always appear in the control region (Control line) if the test has been performed properly.

4. Reagents Provided

Reagents (Store at 4-30°C)	E2810	E2830
Reaction device	10pcs	30pcs
Prefilled sample tube with screw cap (5ml)	10pcs	30pcs
Disposable pipettes	10pcs	30pcs
Prefilled Extraction tube with dropper tip (0.5ml)	10pcs	30pcs
Sterile swab	10pcs	30pcs
Instruction manual	1	1

5. Materials required but not provided

- A grinder sufficient to render sample to particle size of fine instant coffee
- Balance with 0 - 50 g measuring capability
- Microcentrifuge and centrifugal vials
- Vortex mixer and/or Shaker
- S-Flow Reader or 3PRMini along with matching software (for quantitative detection)

6. Storage Instructions

Store kit reagents between 4 and 30°C (39.2 - 86°F). Do not freeze any components provided. Expiry of the kit and reagents is stated on the labels respectively and no quality guarantee is accepted after the expiration date. The expiry of the kit components can only be guaranteed if the components are stored properly as well as if the reagent is not contaminated by the first handling, in case of repeated use of one component. Do not interchange individual reagents between kits of different lot numbers.

7. Safety and Precautions for use

- Use gloves and disinfect the workbench before starting.
- All reagents should be warmed in room temperature before use and covered when not in use. Use a clean disposable pipette for each sample, in order to avoid cross-contamination.
- Clean surfaces, glass vials, mincers and other equipment before and after each sample preparation.
- Do not mix and interchange different samples.
- Do not interchange individual reagents between kits of different lot numbers
- Do not re-use any of the kit components as they are single-use only
- Do not eat or drink in the area where the samples and the kit are stored and handled.

8. Samples Preparation

8.1 Solid Samples

1. The sample must be collected according to established sampling techniques. Grind a representative sample (at least 5 g) to the particle size of fine instant coffee (50% passes through a 20 mesh screen).
2. Weigh out a 0.5 g ground portion of the sample, add it into the prefilled sample tube and vortex it for 1 min. Alternatively, shake vigorously by hand. **The ratio of sample to extraction solvent is 1:10 (w/v).**
3. Two (2) ml of the extract should be centrifuged at high speed for 2 min in reaction caps by using a microcentrifuge. Alternatively, let the sample settle down.
4. Using a disposable pipette, transfer 3 drops from the supernatant to the reaction device and allow the test to develop for 5 minutes.

NOTE 1: The extracted sample should have a pH value of 7.2 - 8.5. If the pH is less than 7.2 or more than 8.5, the pH should be neutralized using NaOH or HCl.

NOTE 2: In case of cloudy, thick samples, that do not allow the mixture to develop, an extra dilution 1:1 with the Extraction Buffer is required before transferring 3 drops to the reaction device. In this case, multiply the final milk protein ppm result x 2.

8.2 Liquid Samples and CIP Solutions

1. Use 0.5 mL of the sample, add it into the prefilled sample tube and vortex it for 1 min. Follow the rest of the procedure as in step **8.1.3**.

NOTE:You can skip using the microcentrifuge, unless your sample is a viscous liquid.

8.3 Surfaces and swab sampling

1. Mark out a swabbing area of approximately 10 x 10 cm.
2. Moisten a swab by dipping into the prefilled extraction tube.
3. Gather the sample with the swab by using a crosshatch technique (**Figure 1.**). Move the swab horizontally, vertically, diagonally while rotating the tip. Repeat this starting from a different angle each time.
4. After the sample collection, place the swab in the extraction tube, rotate the swab forcefully against the side of the tube for 1min. Best results are obtained when the sample is vigorously extracted in the solution. Shake it on a vortex for just a second or you can just tip it with your fingers 3-4 times for better homogenization.
5. Close the extraction tube with the dropper tip. Add 3 drops in the circular window of the reaction device and allow the test to develop for 5 minutes.

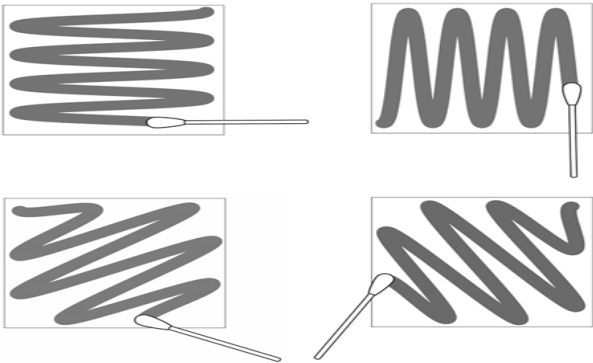


Figure 1.

9. Interpretation of results

Note: For internal procedure purposes three colored lines are present on the result window of the Total Milk Free Test. The colored lines have no effect on the product's performance since they are washed away during the experiment.*

After 5 minutes, the test device can be visually read and interpreted according to the following figure. Observation after 10 minutes lead to inaccurate conclusions.


Negative Result: Two visible colored bands appear at both Hook line (H) and Control (C) line. It indicates that the concentration of milk proteins is zero or below the detection limit of the test.

Positive Result: Three visible colored bands appear at Hook line (H), Control line (C) and Test line (T). Any intensity of Test line indicates the presence of milk proteins into the sample.

High Positive Result: No colored band is visible at Hook line (H) and the band at Test line (T) may be faint or absent. It indicates that the sample contains milk proteins to very high concentrations.*

Invalid Results: No colored band appears at Control line no matter whether it appears at Test line, a Hook line or not.

*In this case you can confirm the high content of your sample by diluting it. Add one drop with the disposable pipette– for solid/liquid samples (or with the dropper tip for surfaces) in a new sample tube with screw cap and shake vigorously. You are expected to get an image as **Positive Result** or **High Positive Result** with the presence of Test line. To quantify your result multiply by the dilution factor x 200.

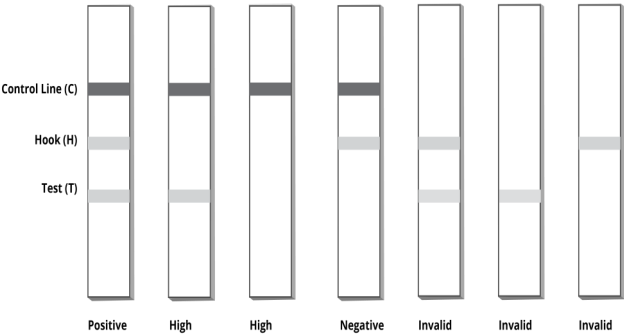

Visual result interpretation index

9.2 Quantitative procedure

After 5 minutes, place the reaction device inside the plastic holder in order to be scanned. Use the appropriate software to quantify results as soon as possible and no later than 1 minute after the end of analysis. The software will use a Lot specific curve to calculate the results. Refer to the Reader's manual for a detailed description of the quantification procedure.

10. Performance Evaluation

10.1 Reference Materials

Several reference materials are being used for the evaluation of each product of ProGnosis Biotech S.A. in the context of Quality Control performed by the Quality Control Department. Please request a validation report, including the results, at exports@prognosis-biotech.com.

11. Assay Claims

- Samples showing negative results may contain milk proteins below the limit of detection of the assay. This Lateral Flow kit does not claim that food is safe for consumption based upon a determination of milk content. Matrix effects may also affect the result of the method.
- The recovery/cross reactivity of the method might be affected when analyzing processed food (e.g. heat treatment, dehydration, etc.), because proteins may be altered or fragmented.
- Food samples that have been heat treated may contain denatured proteins which may not be captured by the antibody. Recovery of these matrices might be reduced.
- The protein content and the protein composition may differ among various species of the same matrix. Therefore, different varieties may produce different results.
- LOD and LOQ in CIP solutions refer to the final rinse water. The presence of cleaning agents and detergents may affect the result of the method.
- Different milk components like whey (β -lactoglobulin), caseinates, etc. are added to food products. These components may affect the milk protein ratio that we used to calibrate this Allergen Free Test. Therefore, if the analyzed food contains whey (β -lactoglobulin) or caseinates then the results may be underestimated.
- Cow's milk contains 3.2% proteins. To express results as ppm milk, multiply the ppm milk proteins result by 31.25 (e.g., 2.5 ppm milk proteins $\times$ 31.25 = 78.1 ppm milk).
- To express results as ppm milk powder, multiply the ppm milk result by 0.1 (e.g., 78.1 ppm milk $\times$ 0.1 = 7.81 ppm milk powder).

12. Method Summary

Total procedure time (after samples and reagents preparation): 5 min.

12.1 Food samples and CIP solutions

Add 0.5gr or 0.5mL of the sample into the prefilled sample Tube



Vortex for 1min or shake by hand

Centrifuge the sample for 2 min, at high speed in a microcentrifuge. Alternatively, let the sample settle down

Transfer 3 drops from the supernatant to the reaction device

Allow test to develop for 5 minutes

Read the results visually or place the device in Reader to be scanned

12.2 Working surfaces

Mark out a swabbing area of approximately 10 x 10 cm



Moisten a swab by dipping into the extraction tube



Gather the sample with the swab by using a crosshatch technique



Place the swab in the prefilled tube to extract the sample



Close the extraction tube with the dropper tip. Add 3 drops in the circular window of the reaction device



Allow test to develop for 5 minutes



Read the results visually or place the device in Reader to be scanned

All immune assays supplied by ProGnosis Biotech S.A., are warranted to meet or exceed our published specification when used under normal conditions in your laboratory. If the product fails during the stated period, a replacement product will be issued.

ProGnosis Biotech S.A. makes no warranty of any kind, either expressed or implied, except that the materials from which its products are made are of standard quality. There is no warranty of merchantability of this product, or of the fitness of the product for any purpose. ProGnosis Biotech S.A. shall not be liable for any damages, including special or consequential damage, or expense arising directly or indirectly from the use of this product. This method is considered to be a screening method, before a legal action, samples detected as positives must be confirmed with a confirmation method. This product is meant to be used only For Research or Manufacturing use and by qualified technicians.

E2810-E2830 Manual_Allergen Free Test_Total Milk_v4_en



LATERAL FLOW TEST KIT

for the detection of milk proteins in food, cip solutions and working surfaces

ProGnosis Biotech S.A. is ISO 9001:2015 certified by TÜV Hellas (TÜV NORD).

Use only the current version of Product Data Sheet enclosed with the kit.

Total Milk Free test, E2810/E2830, is a lateral flow test that detects milk proteins in food products, cip solutions and working surfaces. The lateral flow kit contains all reagents required for the immunoassay method.

Matrices:

Bakery products, beverages, bread mix, cake mix, chocolate, cooked biscuit, dressings, ice cream, sauces, sausage, sorbets, soup mix

- Sample preparation: extraction /swab sampling
- Test time (incubation time after samples and reagents preparation): 5 min
- Calibration Range: 4-100ppm
: 1-32 μ g/100cm² (for surfaces)
- Shelf life: 12 months
- Storage: 4-30°C

Note: The same extract can, also, be used in all the products of Allergen Free series. If you want to check your sample for more than one allergen at the same time, you can use the same extract to test as many of the products you have available.



www.prognosis-biotech.com
e: exports@prognosis-biotech.com
t: +30 2410 623922
Farsalon 153 | 41335 Larissa, Greece

