


# Lovibond<sup>®</sup> Colour Measurement Tintometer<sup>®</sup> Group



## LC 100 Spectrocolorimeter

Hand-Held Imaging Spectrocolorimeter

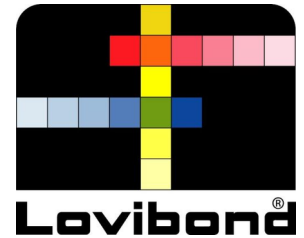

- Handheld and Lightweight Spectrocolorimeter providing Fast, Precise and Accurate Colour Measurement
- Up to 350 Sample Measurements can be Stored and Compared to Reference Standards
- Conformance Reports Available in PDF and Excel Formats
- Precise Positioning on Sample with Preview of Measurement Area

Part Number: 401100

### Handheld and Portable

The LC 100 offers the ultimate in portability; its compact and lightweight design ensures the instrument can be with the user at all times, offering immediate access to accurate colour measurement. Fitting in the palm of your hand the LC 100 ensures easy one-handed operation anywhere on site for a prolonged period of time.

### Accuracy

Despite its lightweight and compact design, the LC 100 does not compromise on accuracy. Using 8 different visible illuminations and 1 UV LED along with the 45°/0 optical geometry and image capture technology, the user can feel confident in the results they are receiving.

### Reporting Results

To make reporting as quick and easy as possible, the LC 100 automatically creates conformance reports in either PDF or Excel formats. The user can then share these with a PC or directly to a compatible printer.

### Storage

The LC 100 can store up to 350 sample measurements, these samples can then be tagged via the LC 100's on-screen keyboard or built in voice recorder. The LC 100 can also store up to 30 reference standards by date and time, ensuring a high level of quality control when comparing results.

### Suitable for Multiple Applications

The versatility of the LC 100 means its suitable for multiple applications and has the ability to report with a display of  $\Delta E$  for all common methods i.e. CIELAB, CMC, CIE 94, 2000, plus a display of measured and delta differences with  $L^*a^*b^*C^*h^\circ$  for standard and sample, graphical  $L^*a^*b^*$  plot.

## Industry

Chemical Industry | Food and Beverage Industry | Industries Others | Oil and Gas | Pharmaceutical Industry

## Application

Chemicals | Coatings | Food and Beverage | Household Products | Medical and Clinical Products | Pharmaceuticals and Cosmetics | Plastics | Sugar Solutions, Syrups and Honey | Textiles

## LC 100 Spectrocolorimeter

The LC 100 is a highly flexible and lightweight spectrocolorimeter. This handheld unit offers the user the ability to capture fast and accurate colour measurement throughout the manufacturing process. In just over a second, the LC 100 will provide a sample reading along with a comprehensive report of how close the colour is to a self-programmed standard.

## Technical Data

<b>Applicable Standard</b>	ASTM E308
<b>Colour Type</b>	Opaque
<b>Display</b>	4.5 cm Colour TFT
<b>Interfaces</b>	USB
<b>Measuring Geometries</b>	45/0 Image Capture
<b>Internal Storage</b>	30 standards/350 sample results
<b>Measurement Time</b>	1.8 seconds
<b>Storage Temperature</b>	-20° to 60°C (-4°F to 104°F)
<b>Short Term Repeatability</b>	Typical 0.10 ΔE 94 on white (D65/10)
<b>Illuminants</b>	CIE Illuminant A CIE Illuminant D65
<b>Portability</b>	Portable

## Delivery Scope

- LC 100 Handheld Spectrocolorimeter
- Carrying Pouch
- UBS Cable for Power Supply and Data Transfer

### Tintometer GmbH

Lovibond® Water Testing  
Schleefstraße 8-12  
44287 Dortmund  
Tel.: +49 (0)231/94510-0  
sales@lovibond.com  
www.lovibond.com  
Germany

### The Tintometer Limited

Lovibond House  
Sun Rise Way  
Amesbury, SP4 7GR  
Tel.: +44 (0)1980 664800  
Fax: +44 (0)1980 625412  
sales@lovibond.uk  
www.lovibond.com  
UK

### Tintometer China

Room 1001, China Life Tower  
16 Chaoyangmenwai Avenue,  
Beijing, 100020  
Customer Care China Tel.: 4009021628  
Tel.: +86 10 85251111 App. 330  
Fax: +86 10 85251001  
chinaoffice@tintometer.com  
www.lovibond.com  
China

### Tintometer South East Asia

Unit B-3-12, BBT One Boulevard,  
Lebu Nilam 2, Bandar Bukit Tinggi,  
Klang, 41200, Selangor D.E  
Tel.: +60 (0)3 3325 2285/6  
Fax: +60 (0)3 3325 2287  
lovibond.asia@tintometer.com  
www.lovibond.com  
Malaysia

### Tintometer Brazil

Caixa Postal: 271  
CEP: 13201-970  
Jundiaí – SP  
Tel.: +55 (11) 3230-6410  
sales@lovibond.us  
www.lovibond.com.br  
Brazil

### Tintometer Inc.

6456 Parkland Drive  
Sarasota, FL 34243  
Tel: 941.756.6410  
Fax: 941.727.9654  
sales@lovibond.us  
www.lovibond.us  
USA

### Tintometer India Pvt. Ltd.

Door No: 7-2-C-14, 2<sup>nd</sup>, 3<sup>rd</sup> & 4<sup>th</sup> Floor  
Sanathnagar Industrial Estate,  
Hyderabad, 500018  
Telangana  
Tel: +91 (0) 40 23883300  
Toll Free: 1 800 599 3891/ 3892  
indiaoffice@lovibond.in  
www.lovibondwater.in  
India

### Tintometer Spain

Postbox: 24047  
08080 Barcelona  
Tel.: +34 661 606 770  
sales@tintometer.es  
www.lovibond.com  
Spain

Technical changes without notice

Printed in Germany

Lovibond® and Tintometer® are Trademarks of the Tintometer Group of Companies