


Lovibond[®] Colour Measurement Tintometer[®] Group



RT520 (diffuse/8) 4 & 8 mm aperture

Flexible Sphere (Diffuse/8°) Spectrophotometer

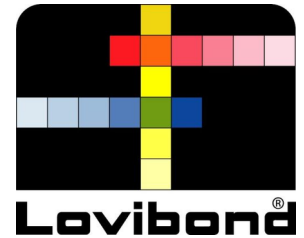

- Switchable 4 mm & 8 mm Measurement Area / 6.5 & 13 mm Target Window
- Fast & Easy Colour Measurement with High-Resolution Colour Screen
- Flexibility with Flip Show Mechanism for Difficult to Reach Areas
- Simple Lock Button Eliminates Mistakes

Part Number: 400520

Improved Profitability

Consistent colour measurements can improve your profitability by avoiding any inconsistencies between customer specification and final production as well as from batch to batch.

Flexible and Configurable to Suit Your Needs

The RT 520 has a large, high-resolution colour screen so you can quickly view the data on-site. The "Project" function allows you to determine a different group of standards for each project and the "Jobs" function gives you step by step instructions, allowing you a programmed sequence of specific steps to assist you and the team in the colour measurement process

Versatility

The flip shoe mechanism adds even more flexibility and makes it easy for you to measure difficult-to-reach areas letting you accurately measure reflective or uneven surfaces or materials including optical brightening agents for a wide variety of products and packaging types. The RT 520 also supports simultaneous SPIN/SPEX output and Correlated Gloss measurements.

Avoid Making Mistakes

A simple lock button ensures that mistakes are not made during quality checks. A remote measurement trigger also allows you to initiate hands-free measurement.

Enhanced Quality Software Control

The RT 520 can be combined with the OnColor[™] software package allowing you to download measurement data, analyse in detail and store indefinitely on a PC.

Industry

Chemical Industry | Food and Beverage Industry | Industries Others | Oil and Gas | Pharmaceutical Industry

Application

Chemicals | Coatings | Food and Beverage | Household Products | Medical and Clinical Products | Pharmaceuticals and Cosmetics | Plastics | Textiles

RT520 (diffuse/8) 4 & 8 mm aperture

The RT 520's large, high-resolution color screen, Graphical Jobs™ interface and configurable menus makes it easy for you to understand colour measurement, bring consistency to your processes and quickly view measurement data.

Technical Data

| | |
|----------------------------|--------------------------------------|
| Applicable Standard | ASTM 313 ASTM D1925 ASTM E308 |
| Colour Type | Transparent Translucent Opaque |
| Spectral Interval | 10 nm - measured; 10nm - output |
| Measurement Range | 0 tot 200 % reflectie |
| Measurement Time | Approx 2 seconds |
| Lamp Life | Approx 500,000 measurements |
| Storage Temperature | -4 to 122°F (-20° to 50°C) |
| Spectral Range | 400 - 700 nm |
| Portability | Benchtop |

Delivery Scope

- Sturdy Box Including
- Instrument
- Calibration box
- Universal power supply
- Instruction manual

Accessories

| Title | Part Number |
|----------------------------------|-------------|
| RT On-Colour Premium Software | 400610 |
| RT300/400/500 Benchtop Stand | 400620 |
| RT300/400/500 Lqd Cell Holder | 400630 |
| RT300/400/500 Powder Cell Holder | 400640 |
| RT V-Block Sample Holder | 400660 |
| RT300/400/500 Pellet Holder | 400680 |

Tintometer GmbH

Lovibond® Water Testing
Schleefstraße 8-12
44287 Dortmund
Tel.: +49 (0)231/94510-0
sales@lovibond.com
www.lovibond.com
Germany

The Tintometer Limited

Lovibond House
Sun Rise Way
Amesbury, SP4 7GR
Tel.: +44 (0)1980 664800
Fax: +44 (0)1980 625412
sales@lovibond.uk
www.lovibond.com
UK

Tintometer China

Room 1001, China Life Tower
16 Chaoyangmenwai Avenue,
Beijing, 100020
Customer Care China Tel.: 4009021628
Tel.: +86 10 85251111 App. 330
Fax: +86 10 85251001
chinaoffice@tintometer.com
www.lovibond.com
China

Tintometer South East Asia

Unit B-3-12, BBT One Boulevard,
Lebuh Nilam 2, Bandar Bukit Tinggi,
Klang, 41200, Selangor D.E
Tel.: +60 (0)3 3325 2285/6
Fax: +60 (0)3 3325 2287
lovibond.asia@tintometer.com
www.lovibond.com
Malaysia

Tintometer Brazil

Caixa Postal: 271
CEP: 13201-970
Jundiaí – SP
Tel.: +55 (11) 3230-6410
sales@lovibond.us
www.lovibond.com.br
Brazil

Tintometer Inc.

6456 Parkland Drive
Sarasota, FL 34243
Tel: 941.756.6410
Fax: 941.727.9654
sales@lovibond.us
www.lovibond.us
USA

Tintometer India Pvt. Ltd.

Door No: 7-2-C-14, 2nd, 3rd & 4th Floor
Sanathnagar Industrial Estate,
Hyderabad, 500018
Telangana
Tel: +91 (0) 40 23883300
Toll Free: 1 800 599 3891/ 3892
indiaoffice@lovibond.in
www.lovibondwater.in
India

Tintometer Spain

Postbox: 24047
08080 Barcelona
Tel.: +34 661 606 770
sales@tintometer.es
www.lovibond.com
Spain