

Haier Biomedical Heating Cabinet

Fast and safe heating without overheating and noise
Economical standard equipment for modern medical facility



Application

Application includes heating application of blood, blood products, liquid medicine, nutrient solution, rinsing liquid and physiological saline. Installations can be found in operating rooms, ICU, emergency rooms, wards, or other related premises.

Produce Features

- Fast heating speed: variable power PTC heating module with fast heating speed and energy saving;
- UV sterilization: UV light sterilization in the cabinet with a variety of optional sterilization methods for effective prevention of possible contamination at the access of the items;
- Multiple overheat protection: electronic control + mechanical temperature control + overheat protection, multiple temperature control to ensure heating safety;
- Environmentally friendly: no compressor, external fan, extremely quiet and low noise with no disturbance to the airflow in the operating room and no effect on temperature, and prevents cross-contamination.

Product Advantages

- PTC ceramic heating modules for fast and efficient heating;
- Standard UV sterilization for effective sterilization; optional sterilization options are available;
- Multiple over heating protections;
- Two-sensor and two-control system
- High temperature and high current protection for heating element

Ergonomics

- Extremely low sound level at 30 db;
- 7-inch LCD touch screen for easier user interaction and information clarity;
- Stainless steel interior and flat external surface are easy to clean;
- Multiple alarms for over-heating protection




Internal Structure4



7 Inch LCD Screen



UV Sterilization Lamp

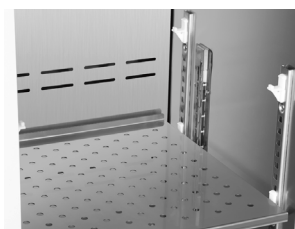

Stainless Steel Partition

Specifications

Model	Voltage (V/Hz)	Power (W)	External Dimension (width*depth*height mm)	Effective Volume	Internal Temperature	Sensor Data	Alarm Mode	Shelf/Storage Basket
HYR- 111	220V/50Hz	130W	597*682*833	111L	Circumstance temperature~50℃	2	Sound + Display	2
HYR- 351	220V/50Hz	360W	597*682*1813	351L	Circumstance temperature~50℃	2	Sound + Display	6

Qingdao Haier Biomedical Co.,Ltd.

Brand Building, Haier Industrial Park, No. 1
Haier Road, Qingdao, 266101, P.R.China
Tel: +86-0532-88935955
Website: www.haiermedical.com

